

# ACE

Veterinary Supplies




Tel: 01491 527027

Email: [sales@acevetsupplies.co.uk](mailto:sales@acevetsupplies.co.uk)

Web: [www.acevetsupplies.co.uk](http://www.acevetsupplies.co.uk)


Issue 8

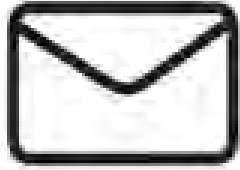

Email: [sales@acevetsupplies.co.uk](mailto:sales@acevetsupplies.co.uk)



Telephone: 01491 527027

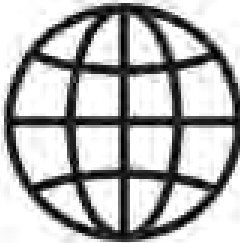

Website: [www.acevetsupplies.co.uk](http://www.acevetsupplies.co.uk)

Find us on Facebook: [@AceVeterinarySupplies](https://www.facebook.com/AceVeterinarySupplies)



Find us on Instagram: [@ace\\_vet\\_supplies](https://www.instagram.com/ace_vet_supplies)



# INDEX

- A** 12 A G Cuffill  
78 Aneroid Sphygmomanometer
- 57 Arterial Catheters
- 3 Avian ET Tubes
- 21 Ayres T Piece
- B** 22 Bain Circuit  
61 Bionector  
63 Biopsy Needles  
55 Bodok Seals
- 11 Bougies
- 32 Breathing Bags
- 20 Breathing Filters
- 17 Broncho-Alveolar Lavage Catheter
- 52 Buster Bag
- C** 86 Capnography Lines  
27 Circle Absorber, SafeAbsorb  
23 Circle Circuits  
8 Cleaning Brushes, Endotracheal  
41 Connectors - Elbows, Luer  
39 Connectors - Straight  
40 Connectors - Straight  
10 Connectors Endotracheal, Low Dead Space  
42 Connectors T & Y  
52 Cradles - Remark  
12 Cuff Pressure Manometer  
62 Curos Disinfecting Port  
55 Cylinder Keys
- D** 72 Defibrillator Pads  
88 D-Lite Connectors  
59 Dog Catheters  
50 Dormasface Electric Blanket  
84 Ear Thermometer
- E** 73 ECG Adapters  
76 ECG Cables - Fukuda Densi  
74 ECG Cables - GE/Datex  
75 ECG Cables - Mindray  
76 ECG Cables -Nihon Kohden  
73 ECG Clips  
73 ECG Gel  
71 ECG Pads  
72 ECG Pads  
65 Edan H100B Veterinary Pulse Oximeter  
15 Eldridge Facemask  
31 EMD Safety Device  
17 Endobronchial Blocker
- 8 Endotracheal Cleaning Brushes  
10 Endotracheal Tube Connectors  
9 Endotracheal Tube Rack
- E** 11 Endotracheal Tube Ties  
8 Endotracheal Tubes Magill
- 4 Endotracheal Tubes PVC Cuffed
- 6 Endotracheal Tubes PVC Cuffed with Suction Lumen
- 6 Endotracheal Tubes PVC Cuffed, K9
- 5 Endotracheal Tubes PVC Cuffed, Portex
- 6 Endotracheal Tubes PVC Cuffed, Reinforced
- 5 Endotracheal Tubes PVC Cuffed, Ventiseal
- 3 Endotracheal Tubes PVC Uncuffed
- 7 Endotracheal Tubes Silicone, Foal
- 7 Endotracheal Tubes Silicone, Large Animal
- 7 Endotracheal Tubes Silicone, Swine
- F** 14 Facemask Diaphragms  
15 Facemask, Eldridge
- 15 Facemask, Oxygen Recovery
- 14 Facemasks, Anaesthetic
- 54 Flowmeters, Oxygen
- 7 Foal Endotracheal Tubes
- 92 FOB Watches
- 48 Foil Bnakets
- 58 Foley Catheters
- 9 Funnel, Replacement, Large Animal ET Tube
- 27 Funnel, Soda Lime
- 86 Gas Sample Lines
- G** 57 Guide Wires
- 50 Heated Table
- H** 18 HMEF's  
19 HME's  
47 Hover Heat
- 51 HUG-U-VAC
- 54 Humidifier Bottles
- I** 80 IBP Transducer Cables  
81 IBP Transducer Cables  
82 IBP Transducer Kits  
89 ICON Oxygen Analyser
- 55 ICU Cage
- 38 Induction chambers
- 11 Intubation Stylets
- 27 IPPV Button Valve
- 60 IV Pressure Bags
- K** 56 Jugular Catheters  
59 KatKath Catheter
- 46 Kennel Connector
- L** 96 Labels  
21 Lack Circuit
- 34 Laryngoscopes - Conventional  
34 Laryngoscopes - Disposable  
33 Laryngoscopes - Miller

<b>L</b>	35	Laryngoscopes - MRI
	33	Laryngoscopes - Veterinary Fibre Optic
	59	Little Herbert
	8	LMA's
	10	Low Dead Space Connectors
<b>M</b>	60	Mantis IV Hanger
	29	Mask, Exotic
	22	Mini Lack Manifold, Darvall
	22	Mini Lack Tubing
	83	Minitemp
	95	MRI Products
	13	Mucosal Device
<b>N</b>	53	Nasal Cannulas
	53	Nasal Oxygen Catheter
	13	Nasal Tube, Anaesthetic
	44	Nervice Blocks
	79	NIBP Cables
	78	NIPB Accessoris
	77	NIPB Cuffs
<b>O</b>	55	O Rings
	91	Oesophageal Stethoscopes
	12	Optilube Tube Lubrication
	89	Oxygen Analyser (ICON)
	90	Oxygen Cell
	54	Oxygen Flowmeters
	53	Oxygen Link Tubing
	53	Oxygen Link Tubing Connector
	15	Oxygen Recovery Facemask
	54	Oxygen Splitters
	54	Oxygen Tube Nipple
<b>P</b>	64	PalmSat VET Pulse Oximeter
	52	Patient Positioning - Buster Bag
	51	Patient Positioning - HUG-U-VAC
	52	Patient Transfer System - Schure Slide
	45	Patient Warming - Mistral Warm Air Heater
	46	Patient Warming Blankets
	47	Patient Warming Blankets - Reusable
	30	PEEP Valves
	30	Pressure Manometer, In circuit
	93	Purair 80 AGSS Systems
	73	Quadrode MRI ECG Pads
<b>R</b>	52	Redmark Cradles
	48	Reflect-A-Therm Table Cover
	43	Regional Anaesthesia
	55	Regulators
	9	Repair Adhesive
	9	Repair Cuffs
	9	Repair Inflation Lines
	16	Resuscitators
<b>S</b>	92	Scavenging APL Valve
	92	Scavenging Charcoal Canisters
	92	Scavenging Connectors

<b>S</b>	93	Scavenging Receivers
	52	Schure Slide Patient Transfer System
	62	Sedation Kit
	12	Silko Spray
	48	Silver Heat Products
	28	Soda Lime, Spherasorb, IS Can & Drum
	66	SpO2 Sensor Clips
	70	SpO2 Sensors - Fukuda Denshi
	69	SpO2 Sensors - GE/Datex
	70	SpO2 Sensors - Masimo
	67	SpO2 Sensors - Mindray
	68	SpO2 Sensors - Nellcor
	68	SpO2 Sensors - Nonin
	66	SpO2 Sensors - Surgivet
	70	SpO2 Sensors - Viamed
	78	Sphygmomanometer
	88	Spirometry Sets
	90	Suction Catheters
	91	Suction Controller
	91	Suction Filter, Eschmann
	90	Suction Jars
	91	Suction Pump, Manual
	90	Suction Tubing
	7	Swine Endotracheal Tubes
	96	Syringe Labels
<b>T</b>	21	T Piece, Vetronics Tiny
	85	Temperature Probes
	49	Thermal Mattresses
	49	Thermal Mattresses - MRI
	84	Thermometer, Ear
	84	Thermometers
	62	Three Way Connector, BD
	36	Tracheostomy
	37	Tracheostomy
	54	Tube Nipple
	9	Tube Rack
	11	Tube Ties, Endotracheal
	58	Urinary Catheters
<b>V</b>	62	V Kit Sedation Kit
	94	Vaporiser Easy Fillers
	94	Vaporiser Key Fillers
	94	Vaporiser 'O' Rings
	94	Vaporiser Pour Fillers
	16	Venturi Valves
	87	Water Traps
<b>Y</b>	25	Y Circuit Smooth Bore
	24	Y Circuit with Water Trap
	24	Y Circuit, Compact (3 metres)
	26	Y Circuit, Darvall
	91	Yankauer Handles
	29	Zero Dead Space Mask

## Avian/Cole Tubes, Silicone

- Sterile
- Radiopaque
- Accepts 15mm Adapter

Also referred to as the 'Cole' Endotracheal Tube, with lower profile on the distal end to aid proper placement and larger proximal end helps form a tighter seal. Stepped down 'Cole' style wall. Complete with standard 15mm connection.

"We find that these tubes meet our needs perfectly, providing both flexibility and rigidity"

Matthew Rendle RVN



Code	Description
ATS010	1.0mm x 12cm with Guide Wire
ATS015	1.5mm x 12cm with Guide Wire
ATS2013	2.0mm ID x 5.0mm x 13cm
ATS2513	2.5mm ID x 6.3mm x 13cm
ATS3015	3.0mm ID x 8.0mm x 15cm
ATS3515	3.5mm ID x 8.6mm x 15cm
ATS4015	4.0mm ID x 9.6mm x 15cm
ATS5018	5.0mm ID x 12.6mm x 18cm
ATS6018	6.0mm ID x 15.1mm x 18cm

## Endotracheal Tubes Uncuffed, PVC

Available individually or in boxes of 10

- Soft Medical Grade PVC
- Transparent for easy observation
- Clearly Marked Guidelines & Depth Markings
- Radiopaque positional guideline
- Murphy Eye

These ET Tubes are available up to 10mm upon request. Maybe subject to minimum quantity of 10.



Single Code	Size	Box Quantity	Box Code
ETUP020	2.0mm	10	ETUP020/10
ETUP025	2.5mm	10	ETUP025/10
ETUP030	3.0mm	10	ETUP030/10
ETUP035	3.5mm	10	ETUP035/10
ETUP040	4.0mm	10	ETUP040/10
ETUP045	4.5mm	10	ETUP045/10
ETUP050	5.0mm	10	ETUP050/10
ETUP055	5.5mm	10	ETUP055/10
ETUP060	6.0mm	10	ETUP060/10
ETUP065	6.5mm	10	ETUP065/10
ETUP070	7.0mm	10	ETUP070/10


01491 527027

## Endotracheal Tubes Cuffed, PVC

Available in singles or in boxes of 10

Code	Size	Box Quantity	Box Code
ETCP025	2.5mm	10	ETCP025/10
ETCP030	3.0mm	10	ETCP030/10
ETCP035	3.5mm	10	ETCP035/10
ETCP040	4.0mm	10	ETCP040/10
ETCP045	4.5mm	10	ETCP045/10
ETCP050	5.0mm	10	ETCP050/10
ETCP055	5.5mm	10	ETCP055/10
ETCP060	6.0mm	10	ETCP060/10
ETCP065	6.5mm	10	ETCP065/10
ETCP070	7.0mm	10	ETCP070/10
ETCP075	7.5mm	10	ETCP075/10
ETCP080	8.0mm	10	ETCP080/10
ETCP085	8.5mm	10	ETCP085/10
ETCP090	9.0mm	10	ETCP090/10
ETCP095	9.5mm	10	ETCP095/10
ETCP100	10.0mm	10	ETCP100/10
ETCP10533	10.5mm x 33cm	10	ETCP10533/10
ETCP1133	11.0mm x 33cm	10	ETCP1133/10

- Soft Medical Grade PVC
- Transparent for easy observation
- Clearly Marked Guidelines & Depth Markings
- Low Inflation Line for Easy Length Reduction
- Radiopaque positional guideline
- Murphy Eye
- Soft Cuff Wall only 8 Microns Thickness

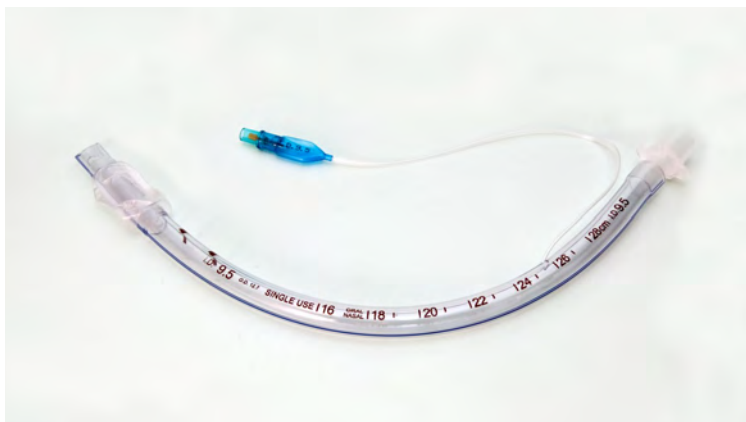

01491 527027


## Endotracheal Tubes Ventiseal™ Cuffed, PVC

Available individually or in boxes of 10

Code	Size	Box Quantity	Box Code
ETCPV025	2.5mm	10	ETCPV025/10
ETCPV030	3.0mm	10	ETCPV030/10
ETCPV035	3.5mm	10	ETCPV035/10
ETCPV040	4.0mm	10	ETCPV040/10
ETCPV045	4.5mm	10	ETCPV045/10
ETCPV050	5.0mm	10	ETCPV050/10
ETCPV055	5.5mm	10	ETCPV055/10
ETCPV060	6.0mm	10	ETCPV060/10
ETCPV065	6.5mm	10	ETCPV065/10
ETCPV070	7.0mm	10	ETCPV070/10
ETCPV075	7.5mm	10	ETCPV075/10
ETCPV080	8.0mm	10	ETCPV080/10
ETCPV085	8.5mm	10	ETCPV085/10
ETCPV090	9.0mm	10	ETCPV090/10
ETCPV095	9.5mm	10	ETCPV095/10
ETCPV100	10.0mm	10	ETCPV100/10
ETCPV105	10.5mm	10	ETCPV105/10
ETCPV110	11.0mm	10	ETCPV110/10

An optimal volume cuff with ideal pressure distribution and effective sealing. The thin ultra-smooth cuff wall reduces abrasive creases and folds, conforms to the tracheal contours and provides an effective seal.

Ensures tube patency for patient safety, whilst softening at body temperature to conform to the patient's respiratory tract.

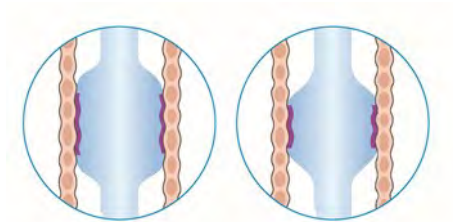

## Endotracheal Tubes Cuffed, Portex

Available in singles or in boxes of 10

Code	Size	Box Quantity	Box Code
ETCPP050	5.0mm	10	ETCPP050/10
ETCPP055	5.5mm	10	ETCPP055/10
ETCPP060	6.0mm	10	ETCPP060/10
ETCPP065	6.5mm	10	ETCPP065/10
ETCPP070	7.0mm	10	ETCPP070/10
ETCPP075	7.5mm	10	ETCPP075/10
ETCPP080	8.0mm	10	ETCPP080/10
ETCPP085	8.5mm	10	ETCPP085/10
ETCPP090	9.0mm	10	ETCPP090/10
ETCPP095	9.5mm	10	ETCPP095/10



Please Note: The Cuffed Portex tubes are not available below 5.0mm.

01491 527027

## Endotracheal Tubes Cuffed, PVC. Larger Sizes

Code	Size
ETCPK912	12.0mm (Curved)
ETCPK913	13.0mm (Curved)



## Endotracheal Tubes Cuffed, with Suction Lumen

- Soft & Flexible Thermosensitive Medical Grade PVC
- Designed with an Evacuation Lumen
- Clearly Marked Guidelines & Depth Markings
- Supplied in Individual Sterile Sachets
- Radiopaque positional guideline
- Murphy Eye

Code	Size	Box Quantity	Box Code
ETCPS060	6.0mm	10	ETCPS060/10
ETCPS065	6.5mm	10	ETCPS065/10
ETCPS070	7.0mm	10	ETCPS070/10
ETCPS075	7.5mm	10	ETCPS075/10
ETCPS080	8.0mm	10	ETCPS080/10
ETCPS085	8.5mm	10	ETCPS085/10
ETCPS090	9.0mm	10	ETCPS090/10



Available individually or  
in boxes of 10

## Endotracheal Tubes Cuffed, Reinforced PVC

Code	Size	Box Quantity	Code	Code	Size	Box Quantity	Code
ETRCP030	3.0mm	5	ETRCP030/5	ETRCP070	7.0mm	5	ETRCP070/5
ETRCP035	3.5mm	5	ETRCP035/5	ETRCP075	7.5mm	5	ETRCP075/5
ETRCP040	4.0mm	5	ETRCP040/5	ETRCP080	8.0mm	5	ETRCP080/5
ETRCP045	4.5mm	5	ETRCP045/5	ETRCP085	8.5mm	5	ETRCP085/5
ETRCP050	5.0mm	5	ETRCP050/5	ETRCP090	9.0mm	5	ETRCP090/5
ETRCP055	5.5mm	5	ETRCP055/5	ETRCP095	9.5mm	10	ETRCP095/10
ETRCP060	6.0mm	5	ETRCP060/5	ETRCP100	10.0mm	10	ETRCP100/10
ETRCP065	6.5mm	5	ETRCP065/5				



These resemble standard endotracheal tubes, but also contain a spiral of metal in the tube wall. This is particularly useful where kinking of the tube may occur perhaps during dental work or eye surgery.

Sizes 3mm — 9.0mm are supplied with an intubation stylet

Available individually or in  
box quantities

**01491 527027**



## Endotracheal Tubes Silicone, Foal

- Tubes up to 12 mm have a 15 mm airway adapter fitting
- Silicone funnel is supplied on the 14 mm tube size
- Repair balloon cuffs available
- Replacement inflation lines available
- Autoclavable
- Radiopaque
- Non-sterile
- 55 cm or 22-inches long



Code	Size ID	Length
FN7	7.0mm	55cm
FN8	8.0mm	55cm
FN9	9.0mm	55cm
FN10	10.0mm	55cm
FN11	11.0mm	55cm
FN12	12.0mm	55cm
FN14	14.0mm	55cm

## Endotracheal Tube, Silicone Swine Extra Long 8mm

Code	Size
ETCSS080	8.0mm ID x 47cm

Especially designed with the co-operation of leading research centres for use in swine. Silicone. Graduated markings with standard adaptor, inflatable cuff, radio-paque, and lateral Murphy eye, extra long length 8mm I.D. x 47cm.

## Endotracheal Tubes Silicone, Large Animal

Code	Size
ETCSE14	14.0mm ID, 19mm OD x 57cm
ETCSE16	16.0mm ID, 22mm OD x 70cm
ETCSE18	18.0mm ID, 24mm OD x 75cm
ETCSE20	20.0mm ID, 26.5mm OD x 90cm
ETCSE22	22.0mm ID, 30mm OD x 90cm

Code	Size
ETCSE24	24.0mm ID, 32mm OD x 90cm
ETCSE26	26.0mm ID, 35mm OD x 100cm
ETCSE28	28.0mm ID, 39mm OD x 100cm
ETCSE30	30.0mm ID, 42mm OD x 100cm

Universal fitting allows connection to most large animal anaesthesia machines and ventilators



01491 527027

## Endotracheal Tubes, Magill Cuffed

Code	Size	Code	Size
<b>ETCM12036</b>	12.0mm x 36cm	<b>ETCM14036</b>	14.0mm x 36cm
<b>ETCM12048</b>	12.0mm x 48cm	<b>ETCM14048</b>	14.0mm x 48cm
<b>ETCM13036</b>	13.0mm x 36cm	<b>ETCM16036</b>	16.0mm x 36cm
<b>ETCM13048</b>	13.0mm x 48cm	<b>ETCM16048</b>	16.0mm x 48cm



Other sizes and uncuffed, are available on request.

For ET Tube Connectors see Page 10

## LMA's



Code	Size	Code	Size
<b>LMA010</b>	Size 1	<b>LMA030</b>	Size 3
<b>LMA015</b>	Size 1.5	<b>LMA040</b>	Size 4
<b>LMA020</b>	Size 2	<b>LMA050</b>	Size 5
<b>LMA025</b>	Size 2.5		

Silicone & Reinforced LMA's are available on request.

Please ask for details and prices.

## Endotracheal Tube Cleaning Brushes

### Large Equine

Code	Size
<b>ETCBEQLA</b>	Large Equine



### Small Animal

Code	Size
<b>ETCBLAR</b>	Large
<b>ETCBMED</b>	Medium
<b>ETCBSML</b>	Small



01491 527027

## Silicone Endotracheal Tube Repair

### Repair Cuffs

Code	Size	Code	Size
RC160	16mm	RC240	24mm
RC180	18mm	RC260	26mm
RC200	20mm	RC280	28mm
RC220	22mm	RC300	30mm



### Repair Cuff Adhesive

Code	Description
RCADH	Repair Cuff Adhesive 82.8ml (85gm)

### Repair Inflation Lines

Code	Size
RILL	For ET Tube sizes 16mm—30mm
RILS	For ET Tube sizes 3mm—14mm



### Replacement Funnel

Code	Size
ETFC32364B1	Large Animal Replacement Blue Funnel Connector



### ET Tube Rack

Code	Description
ETTR	Endotracheal Tube Rack



Acrylic rack for mounting 3 x 3.0mm—14mm ET Tubes Horizontal to aid drainage. 66cm long.

01491 527027

## Endotracheal Tube Connectors

### Small Animal

Code	Size	Code	Size	Code	Size
ETTC025	2.5mm	ETTC055	5.5mm	ETTC085	8.5mm
ETTC030	3.0mm	ETTC060	6.0mm	ETTC090	9.0mm
ETTC035	3.5mm	ETTC065	6.5mm	ETTC095	9.5mm
ETTC040	4.0mm	ETTC070	7.0mm	ETTC100	10.0mm
ETTC045	4.5mm	ETTC075	7.5mm	ETTC105	10.5mm
ETTC050	5.0mm	ETTC080	8.0mm	ETTC110	11.0mm



### Large Animal

Code	Size	Code	Size
ETTC120	12.0mm Metal	ETTC120P	12.0mm Plastic
ETTC130	13.0mm Metal	ETTC130P	13.0mm Plastic
ETTC140	14.0mm Metal	ETTC140P	14.0mm Plastic
ETTC160	16.0mm Metal		



### Low Dead Space

Code	Size
ETTCLD025	2.5mm Low Dead Space
ETTCLD030	3.0mm Low Dead Space
ETTCLD035	3.5mm Low Dead Space
ETTCLD040	4.0mm Low Dead Space
ETTCLD050	5.0mm Low Dead Space



### Connectors with Side Port

Code	Size
SSA2215	Side Stream Adapter 22/15mm Female with Sample Port



### Large Equine

Code	Size	Code	Size
ETTC180	18.0mm	ETTC300	30.0mm
ETTC200	20.0mm	ETTC350	35.0mm
ETTC250	25.0mm	ETTC400	40.0mm



These connectors allow the attachment of the large animal Endotracheal Tube to a small animal anaesthetic circuit.

01491 527027

## ET Tube Ties


Code	Description
<b>ETTT50</b>	ET Tube Ties. Pack of 50



## Intubation Aids

### Stylets

Code	Size	Box Qty	Box Code
<b>IS6</b>	6Fr, 26cm	20	<b>IS6/20</b>
<b>IS10</b>	10Fr, 33cm	20	<b>IS10/20</b>
<b>IS14</b>	14Fr, 33cm	20	<b>IS14/20</b>
<b>IS15</b>	15Fr, 70cm	10	<b>IS15/20</b>

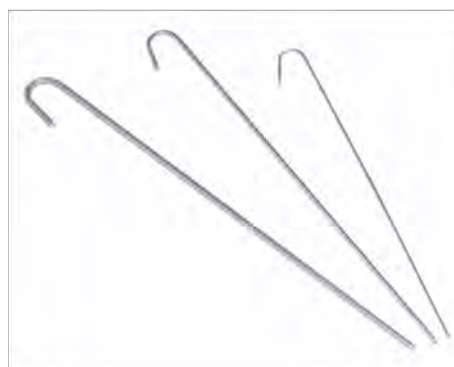

Available individually or in quantity boxes

Malleable aluminium with a low friction plastic sheath to aid insertion of ET Tubes.

Fully rounded atraumatic tip to minimum risk of patient trauma

Sterile and individually packed.

## Bougies

For directional control during routine or difficult endotracheal intubation when the laryngeal inlet cannot be completely seen.

### Vented

Code	Size
<b>IBV10</b>	10Fr, 70cm
<b>IBV15</b>	15Fr, 70cm

### Standard

Code	Size
<b>IB06</b>	6Fr, 53.3 cm
<b>IB10</b>	10Fr, 70cm
<b>IB15</b>	15Fr, 70cm



Code	Size
<b>IBH800</b>	Bougie Holder

Holds 4 x 700/800mm Bougies. Supplied Empty

01491 527027

## Endotracheal Tube Lubrication

### Optilube Lubricating Jelly



Code	Description
OPTL5150	150 x 5G Sachets

### Silko Spray 500ml

Code	Description
SILSP500	Rusch Silko Spray 500ml



## Cuff Pressure Management

### AG Cuffill

Code	Description
AGCCP	AG Cuffill Device



Allows measuring and reading the exact cuff pressure without losing pressure during connection. Easy to use as the operation of the device is similar to the existing practice of filling the cuff with a syringe. Provides indication of cuff pressure & volume. Can be used approximately 100 times.

### Cuff Pressure Manometer

Code	Description
CICPM	Cuff Pressure Manometer
CICPM-SA2	Replacement 1 metre extension tube

Manual cuff pressure gauge for inflating and controlling tracheal tubes with high volume low pressure cuffs. For continued cuff pressure monitoring, the device can be connected to the tracheal tube by means of a 1m extension tube.



01491 527027

## Surgivet Post Anaesthetic Nasal Tube

Code	Description
<b>ENAT14</b>	Nasal Tube 14.0mm ID x 19mm OD x 40cm

- A proximal silicone funnel fitting to avoid tube migration
- Autoclavable
- Non-sterile
- Radiopaque



## Intersurgical Intranasal Mucosal Device



Medication can be quickly and painlessly delivered with the DART Intranasal mucosal atomisation device. Atomised medication delivery results in therapeutic onset delivery similar to intravenous delivery. There is no specialised training required for a system wide adoption. The delivery method is painless for the patient, helping them remain calm in a high stress situation.



Code	Description
<b>DART140</b>	DART™ intranasal atomisation device with 3ml syringe and vial adapter.
<b>DART300</b>	DART300™ intranasal atomisation device without syringe

## Laryngo-tracheal Mucosal Atomisation Device

Medication can be quickly and effectively delivered for difficult airways with the DART-Reach. The device has a flexible stylet, allowing for a customised approach and can be used in a variety of procedures, for example awake intubations, transoesophageal echocardiograms, or difficult airways.



Code	Description
<b>REACH700</b>	DART-Reach™ Atomisation Device, 218mm
<b>REACH720</b>	DART-Reach™ Atomisation Device, 117mm



01491 527027

## Anaesthetic Facemasks

Code	Size
<b>AFMXXS50</b>	XX Small 50mm dia
<b>AFMXS60</b>	X Small 60mm dia
<b>AFMS65</b>	Small 65mm dia
<b>AFMM90</b>	Medium 90mm dia
<b>AFML110</b>	Large 110mm dia
<b>AFMXL140</b>	X Large 140mm dia

Made of clear plastic, with black rubber diaphragm. No latex. Supplied with elastic neck strap (except 50mm & 60mm) for better fixation of the mask. 15 mm connection piece fits inner patient hose.



The connectors below will fit your anaesthetic facemask if the one you have gets lost or broken...

**INT1943** Straight connector  
15M - 22M



**INT1997** Fixed elbow 15M-22M/15F



## Diaphragms

Replacement Diaphragms for Anaesthetic Facemasks

Code	Size
<b>DIAXXS50</b>	XX Small 50mm dia
<b>DIAXS60</b>	X Small 60mm dia
<b>DIAS65</b>	Small 65mm dia
<b>DIAM90</b>	Medium 90mm dia
<b>DIAL110</b>	Large 110mm dia
<b>DIAXL140</b>	X Large 140mm dia



01491 527027



## Eldridge Anaesthetic Facemasks

The Eldridge Anaesthesia Mask addresses a common problem. A standard round mask wants to roll away from patient on a flat table. Towels are often used to brace the mask into position. The Eldridge mask is flat on one side to prevent this "rolling" and safely keep the mask in position. Patent pending.

- Used to induce gas anaesthesia and sometimes maintain anaesthesia in patients
- Provides oxygen therapy in patients with severe respiratory distress



Code	Size
<b>ELDFM-SML</b>	Small, OD 3.5", Depth, 3", Opening 1"
<b>ELDFM-LAR</b>	Large, 5" OD 5", Depth, 4.75", Opening 1.75"
<b>ELDFM-XLA</b>	Extra Large, OD 6", Depth, 5.75", Opening 3"
<b>ELDFM-SET3</b>	Set of 3 Small, Large, X Large

## Eldridge Diaphragms

Code	Size
<b>ELDIA-SML</b>	Small, OD 3.5", Opening 1"
<b>ELDIA-LAR</b>	Large, 5" OD 5", Opening 1.75"
<b>ELDIA-XLA</b>	Extra Large, OD 6", Opening 3"

## Oxygen Recovery Masks

- Facilitate post-surgery, post-trauma and post smoke inhalation patient recovery
- Mask body is manufactured from durable polycarbonate
- Diaphragm is made from tear-resistant thermoplastic material and is inert with polycarbonate for longer life
- Mask features dual vents, not valves, to ensure unrestricted inhalation and exhalation to prevent unforeseen reduction or increase in oxygen flow
- Includes a rubber-mounted 22mm oxygen adapter and 7ft of oxygen tubing
- Designed to be used with oxygen tanks
- Optional hook and strap sets available to secure the mask
- Used by fire departments around the world
- Easily cleaned in 1:20 solution of bleach and water.
- 2 year warranty against manufacturer defects

Code	Size
<b>ORM-SML</b>	Oxygen Recovery Mask — Small
<b>ORM-LAR</b>	Oxygen Recovery Mask — Large
<b>ORM-FEL</b>	Oxygen Recovery Mask — Feline
<b>ORM-PUP</b>	Oxygen Recovery Mask —Puppy/Kitten
<b>ORM-RAP</b>	Oxygen Recovery Mask - Raptor



01491 527027

## Venturi Valves

Code	Size
<b>OX146</b>	Venturi Valve Kit comprising of 24%, 28%, 31%, 35%, 40% and 60%



## Resuscitators

### Semi-Disposable

Code	Size
<b>RES0280</b>	280ml
<b>RES0600</b>	600ml
<b>RES1000</b>	1000ml
<b>RES1650</b>	1650ml



### Autoclavable

The Reusable Resus Bag is manufactured from high quality materials that can be autoclaved up to 50 cycles.

Code	Size
<b>RES0280A</b>	Small PLV 40cm H <sub>2</sub> O
<b>RES0500A</b>	Medium PLV 40cm H <sub>2</sub> O
<b>RES1500A</b>	Large PLV 60cm H <sub>2</sub> O



## Foal/Calf

Code	Description
<b>RES-FOAL</b>	Foal/Calf Resuscitator

- Performance-tested Resuscitator
- Aspirator Mask
- Resuscitator Mask
- Carrying Case
- Instructions
- 2 year warranty



01491 527027

## Broncho-Alveolar Lavage Catheter (Silicone)

- The distal balloon can be inflated to prevent retrograde fluid flow
- The proximal two-way fitting allows simultaneous connection for both irrigation and aspiration
- Repair balloon cuffs and silicone adhesive are available
- Replacement inflation lines are available
- Autoclavable
- Sterile



Code	Tube ID	Tube OD	Needle Length
<b>BAL240</b>	3.0mm	11mm	240cm
<b>BAL300</b>	3.0mm	11mm	300cm

## Endobronchial Blocker

- Can be inserted down a tracheal tube after tracheal intubation to block the right or left main bronchus for a controlled one lung ventilation. The lung tissue distal to the obstruction will collapse, allowing the surgeon's view and access to the relevant structures within the thoracic cavity
- Achieves lung separation and one lung ventilation as an alternative to double lumen tubes
- Manufactured from transparent medical grade PVC allowing clinicians to visualise any blockage quickly
- Atraumatic tip with double offset eye to provide relief if the distal hole becomes occluded for easy insertion
- Smooth surface for easy insertion
- Supplied sterile in easy access sachets



Code	Description
<b>EBO7FR</b>	Endobronchial Blocker 07FR
<b>EBO9FR</b>	Endobronchial Blocker 09FR

01491 527027

## HMEFs

Heat and Moisture Exchanging Filters (HMEFs) combines the efficiency of dedicated breathing filters with optimum moisture return



### Intersurgical Clear-Therm™

Available individually or in box quantities



Single Code	Box Qty	Box Code
HME1441	20	HME1441/20
HME1831	40	HME1831/40
HME1541	150	HME1541/150

#### HME1441 (Micro)

Dead Space 11ml  
 Capnograph Port  
 Connectors 15M/15F  
 Box Quantity 20

#### HME1831 (Midi)

Dead Space 26ml  
 Capnograph Port  
 Connectors 22M/15M-22M/15F  
 Box Quantity 40

#### HME1541

Dead Space 59ml  
 Capnograph Port  
 Connectors 22M/15M-22M/15F  
 Box Quantity 150

## GE Healthcare HMEF with Luer Port

Available individually in box quantities



Code	Box Qty	Box Code
HME2106570-010	50	HME2106570-010/50
HME2106570-009	50	HME2106570-009/50

#### HME2106570-009 (Medium)

Dead Space 25ml  
 Capnograph Port  
 Connectors 22M/15M-22M/15F  
 Box Quantity 50

#### HME2106570-010 (Large)

Dead Space 66ml  
 Capnograph Port  
 Connectors 22M/15M-22M/15F  
 Box Quantity 50



01491 527027

## HMEs

Heat and moisture exchangers (HMEs) designed to replicate the functions of the body's upper airway by conserving expired heat and moisture and returning these to the patient during inhalation.

### Intersurgical Hydro-Therm™



#### HME1442

Dead Space 2.2ml  
Connectors 15M/15F  
Box Quantity 30  
Can be used with ET Tube connectors with Capnograph Connector



#### HME1560

Dead Space 59ml  
Capnograph Port  
Connectors 15M-22M/15F  
Box Quantity 150

Available individually or in box quantities.



#### HME1855

Dead Space 20ml  
Capnograph Port  
Connectors 15M-22M/15F  
Box Quantity 20

Code	Box Qty	Box Code
HME1442	30	HME1442/30
HME1855	20	HME1855/20
HME1560	150	HME1560/150

### Hydro-Trach T HME Filter

A heat, moisture exchanger designed for use on tracheostomised patients. The Hydro-Trach is an ideal product for prolonged use with spontaneously breathing patients.

Available individually or in box quantities

Code	Box Qty	Box Code
INT1875	25	INT1875/25

#### INT1875

Compressible Volume 19ml  
Weight 8g  
Box Quantity 25



01491 527027

## Breathing Filters



### Intersurgical Inter-Guard™

The Inter-Guard™ range of sterile breathing filters is designed for use in breathing systems in the operating theatre for the protection of the patient, breathing system and equipment.



#### BF1344

Dead Space 42ml  
 Capnograph Port  
 Connectors 22F/15M-22M/15F

### Intersurgical Filta-Guard™

The high-efficiency Filta-Guard™ is a dedicated breathing filter designed for use in breathing systems in anaesthesia for the protection of the patient, hospital personnel and the equipment from potential microbial contamination.



#### BF1944

Dead Space 67ml  
 Capnograph Port  
 Connectors 22F-22M/15F  
 Box Quantity 70

### Intersurgical Air-Guard

#### For use in respiratory systems and oxygen concentrators

Air-Guard™ is a hydrophobic pleated membrane filter for the protection of oxygen concentrators and other respiratory equipment.

Providing an excellent level (>99.999%) of protection against bacterial and viral challenges, the Air-Guard's product performance has been tested, validated and verified at independent microbiological laboratories.

The Air-Guard pleated mechanical filter is validated for 24 hours use but this can be extended when used in non-humidified gas flow.

Please note: this product is contraindicated for use at the patient connection end of the breathing system.



#### BF1790

Volume 120ml  
 Connectors 22F-22M/15F  
 Box Quantity 50

Single Code	Box Qty	Box Code
BF1344	50	BF1344/50
BF1944	70	BF1944/70
BF1790	50	BF1790/50

01491 527027

## Ayres T Piece



Available individually or in box quantities

Single Code	Description	Box Qty	Box Code
INT2122	Ayres T Piece with APL Valve & 0.5 litre bag	15	INT2122/15
INT2811	1.0 litre breathing bag with 15mm Neck to Fit Ayres T piece	10	INT2811/10
INT2805	Replacement 0.5 litre breathing bag	10	INT2805/10
INT1661	Replacement End Connector - Intersurgical Straight Connector 15M/11.5-10M	35	INT1661/35



To extend your Ayres T Piece you will need....

Code	Description
INT1576	15mm Corrugated Tubing cut to the required length. Cuffed at 400mm.
INT1962	Connector 15F—22M (for patient limb)
INT1945	Connector 15M-15M (for bag limb)

INT1945



INT1962




## Vetric Tiny T Piece



Code	Description
VTP	Vetric Ultra-low dead-space T-piece circuit for use with patients 500g or less

## Lack Circuit (Mapleson A)

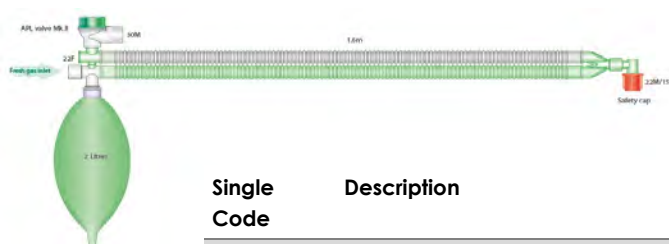

Available individually or in boxes of 10


Single Code	Description	Box Qty	Box Code
CIRPL1622I	Parallel Lack System 1.8M with APL and 2 litre bag 22mm tubing — Intersurgical	10	CIRPL1622I/10

01491 527027

## Bain Circuit (Mapleson D)

Available individually or in box quantities

The Bain circuit is a modification of the Mapleson D system. It is a co-axial system in which the fresh gas flows through a narrow inner tube within the outer corrugated tubing.

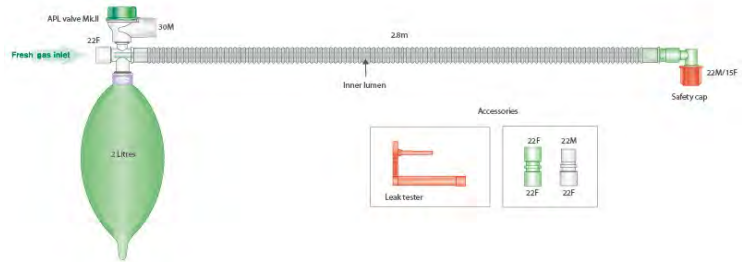

Single Code	Description	Box Qty	Box Code
<b>BAIN2115</b>	Bain 22mm x 1.6m (Intersurgical)	30	<b>BAIN2115/30</b>
<b>BAIN28</b>	Bain 22mm x 2.8m (Intersurgical)	10	<b>BAIN28/10</b>
<b>BAIN54</b>	Bain 22mm x 5.4m (Intersurgical)	5	<b>BAIN54/5</b>



## Mini Lack Tubing

Code	Description
<b>CIRMLC16</b>	Smooth Bore, Mini Lack Y Circuit 1.6m with Y piece (does not include manifold).
<b>CIRMLC16C</b>	Corrugated, Mini Lack Y Circuit 1.6m with Y piece (does not include manifold).



## Darvall Mini Lack Manifold



Single Code	Description
<b>MLACK</b>	Darvall Lack Manifold



01491 527027




## Corrugated Circle ('Y') Circuits

### 22mm Diameter



Intersurgical Flextube breathing system consisting of 2 limbs 22mm diameter, Y-piece, and safety cap. Also contains a 22M straight connector to convert the machine port from 22F to 22M. These circuits can be made to a length of your choice — 3.2 metres and 6 metres already available but please contact our sales office if you need a different length.

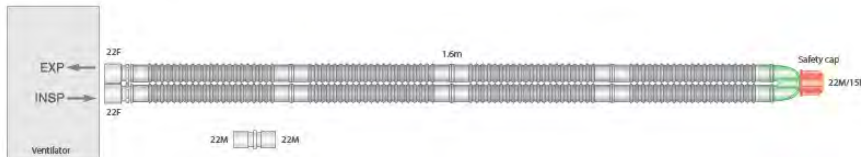


Available individually or in box quantities for certain items

Single-Code	Description	Box Qty	Box Code
<b>CIR2216</b>	Circle Circuit 22mm x 1.6m (2000)	20	<b>CIR2216/20</b>
<b>CIR2224</b>	Circle Circuit 22mm x 2.4m (2001)	15	<b>CIR2225/15</b>
<b>CIR2232</b>	Circle Circuit 22mm x 3.2m		
<b>CIR2260</b>	Circle Circuit 22mm x 6.0m		

### 15mm Diameter

Intersurgical Flextube breathing systems use 15mm tube to reduce circuit volume, especially important when using on the smaller patients. These circuits can be made to a length of your choice — the lengths we currently make are listed below but please contact our sales office if you need a different length.

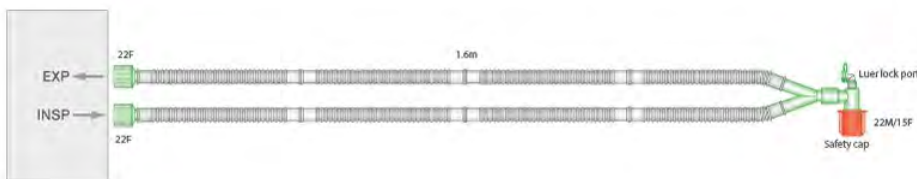


Available individually or in box quantities for certain items

Single-Code	Description	Box Qty	Box Code
<b>CIR1516</b>	Circle Circuit 15mm x 1.6m (2142)	10	<b>CIR1516/10</b>
<b>CIR1520</b>	Circle Circuit 15mm x 2.0m		
<b>CIR1524</b>	Circle Circuit 15mm x 2.4m		
<b>CIR1532</b>	Circle Circuit 15mm x 3.2m		
<b>CIR1560</b>	Circle Circuit 15mm x 6.0m		

01491 527027

## Flextube 'Y' Circuit with Water Trap

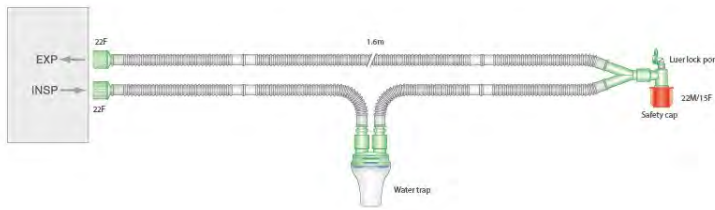

Available individually or in boxes of 10

Single Code	Description	Box Qty	Code
CIR2143	15mm Flextube anaesthetic breathing system with luer elbow and resealable water trap, 1.6m	10	CIR2143/10

## Intersurgical Compact™ 'Y' Circuits

Features and Benefits include:

**Optimum compression** - the unique tube design results in a system which can be compressed to almost one fifth of its extended length.



**Lower resistance to flow** - the Compact tube is designed with a shallow profile resulting in more laminar gas flow and lower resistance to flow.

**Lightweight system** - the Compact system is lighter than other extendable systems ensuring that the 'pull' on the patient's airway device is kept to a minimum.

**Adjustable length** - the Compact system is packed and supplied in the shortest length. This can be simply extended to the length required for clinical use (with the advantage of reduced compressible volume), and compressed for easy disposal. The unique tube has been designed so that it will expand to between four to five times it's packed length, providing the most adjustable system available.

**Easy positioning** - the nature of the Compact system tube allows for it to be set in any position required. This results in multiple positioning possibilities, ideal for head and neck surgery where the system can be positioned away from the patient, and does not require any additional patient connections.

Single Code	Description	Box Qty	Box Code
CIR2158	Circle Circuit 22mm x 3.0m Adjustable	40	CIR2158/40
CIR2163	Circle Circuit 15mm x 3.0m Adjustable	40	CIR2163/40

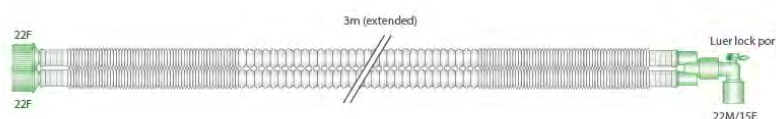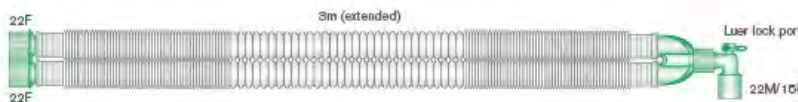

01491 527027

## Make your own Y Circuit to the length you want.....

Should you need a longer Y Circuit we can provide all the parts you need to assemble your own. These include the correct connectors to attach to your Circle System and the Y Connector for the patient end. The 15mm parts list includes a paediatric Y piece to reduce Dead Space. In both sizes the length must be in multiples of 400mm.

To make a 15mm Circuit you will need...



INT1962



INT1902

Code	Description
<b>INT1576</b>	15mm Corrugated Tubing cut to the required length. Cuffed at 400mm. Price per Cuffed Length
<b>INT1962</b>	Connector 15F—22M. 2 required per Circuit
<b>INT1902</b>	15mm Y-piece 15F patient connection - 15M-15M

To make a 22mm Circuit you will need...



INT1961



INT1926

Code	Description
<b>INT1574</b>	22mm Corrugated Tubing cut to the required length. Cuffed at 400m. Price per Cuffed Length
<b>INT1961</b>	Connector 22F—22M. 2 required per Circuit
<b>INT1926</b>	Parallel Y Piece 22MM


### Optional...



Code	Description
<b>INT1978</b>	Safety Cap (Only fits 22mm Y Piece)

## Smooth Bore Circle Circuits & Tubing

Code	Description
<b>CIRSBC15</b>	Smooth Bore Y Circuit 15mm x 1.8m
<b>CIRSBC22</b>	Smooth Bore Y Circuit 22mm x 1.8m
<b>CIRSBC1530</b>	Smooth Bore Y Circuit 15mm x 3.0m
<b>CIRCSBT1530</b>	Smooth Bore Tube 15mm x 3m with 22F end connections
<b>CIR5024</b>	Intersurgical Smooth Bore Y Circuit 15mm x 2.4m

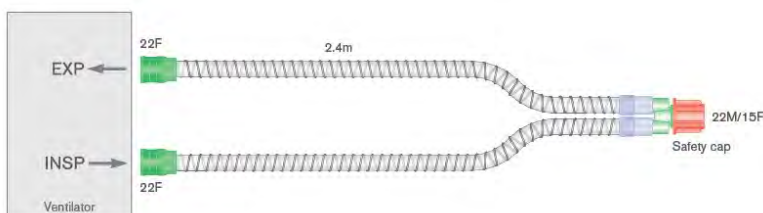

01491 527027

## Darvall Smooth Wall Circle Circuits



**Darvall Smooth Wall Tubing (SWT) Circuits** use low resistance, low volume, small diameter tubes. This new concept in circuit design offers HUGE efficiency advantages for the patient as well as for your practice.

### Features & Benefits:

- Easy to use
- Fits all standard anesthetic gas machines
- Smooth Wall Tubing (SWT) Circuits reduce circuit resistance
- Reduced resistance means they can be used for smaller patients
- Smaller diameter tubing reduces circuit volume by up to 75%
- Reduced circuit volume improves response and efficiency
- SWT Circuits are kink-resistant
- Easy to clean and dries thoroughly and quickly
- Two sizes: Small (12 mm) for animals up to 20kg, and Large (15 mm) for animals over 20kg



Code	Description
------	-------------

<b>CIRD12</b>	12mm x 1.6m, 0.5 & 1 litre bag. For animals 3-20kg
---------------	--

<b>CIRD15</b>	15mm x 1.6m, 3 litre bag. For animals over 20kg
---------------	---

## Corrugated Tubing

Available per 400mm length or in box of 50 metres.

Code per 400mm	Description	Box Qty	Code
<b>INT1574</b>	22mm Corrugated Tubing cut to the required length. Cuffed at 400m. Price per Cuffed Length	50M (125 cuffed lengths)	<b>INT1574-50M</b>
<b>INT1576</b>	15mm Corrugated Tubing cut to the required length. Cuffed at 400mm. Price per Cuffed Length	50M (125 cuffed lengths)	<b>INT1576-50M</b>



## Catheter Mounts

The Superset range of catheter mounts have a concertina design to allow for positioning control. This allows the catheter mount to be set in any required position.

Code	Description	Box Qty	Box Code
<b>INT3514</b>	22/F—22M/15F with Elbow and Luer Lock Port	75	<b>INT3514/75</b>
<b>INT3512</b>	22/F—22M/15F Swivel	75	<b>INT3512/75</b>
<b>INT3505</b>	22/F—22M/15F Swivel with Flip Top Cap	50	<b>INT3505/50</b>
<b>INT3535</b>	15/M - 7.6 port -15/F 46-100mm Catheter Mount paediatric/neonatal	15	<b>INT3535/15</b>



Available individually or in box quantities.

**01491 527027**

## Circle Absorbers

Code	Description
<b>KABAU</b>	Darvall SafeAbsorb (for use with Metal T Piece)
<b>KABTP</b>	22mm T Piece
<b>KABLE</b>	Long Extension for KAB Absorber 22F/22M
<b>SDLFUN</b>	Soda Lime Funnel
<b>KABBM</b>	KAB Bag Mount Elbow (spare)



**KABTP**  
Metal T Piece



**KABBM**  
Spare Bag Mount Elbow 22M-22F



**KABLE**  
Extension for KAB Absorber where, for instance, there needs to be clearance of the vaporiser.



**SDLFUN**  
Powder Funnel to aid the filling of r the KAB Absorber

## I.P.P.V. Button Valve

For use with any breathing system with a 30mm scavenge valve. Just press the button to occlude flow, system instantly re-opens when released.

Code	Description
<b>IPPV-BV</b>	I.P.P.V. Button Valve



01491 527027

## Spherasorb™ Soda Lime

Spherasorb's chemical formulation has been developed specifically to address the potential problems of use within the medical environment.

LoFloSorb eliminates the risk associated with reactions with volatile anaesthetics.



Available individually or in boxes of 2

Code	Description	Box Qty	Box Code
<b>SDLPW45</b>	Spherasorb Pink to White 4.5kg	2	<b>SDLPW45/2</b>

<b>SDLWV45</b>	Spherasorb White to Violet 4.5kg	2	<b>SDLWV45/2</b>
----------------	-------------------------------------	---	------------------

Code	Description	Box Qty	Box Code
<b>SDLLFS45</b>	LoFloSorb Green to Violet 4.5kg	2	<b>SDLLFS45/2</b>

## IS Can™

IS Can, Spherasorb™ disposable CO2 absorber, pink to white or white to violet colour change or Loflosorb™ green to violet colour change. For use with GE® anaesthetic machines: Aisys®, Aespire®, Avance® and ADU®



Code	Description
<b>SDLPW2197</b>	Pink to White Spherasorb® Cannister. Box of 6.

<b>SDLWV2196</b>	White to Violet Spherasorb® Cannister. Box of 6.
------------------	--

<b>SDLGV2198</b>	Green to Violet Loflosorb® Cannister. Box of 6.
------------------	---

## The Drum™



Code	Description
<b>SDLPW2187</b>	Pink to White Spherasorb® 1kg Drum. Box of 10.

<b>SDLWV2186</b>	White to Violet Spherasorb® 1kg Drum. Box of 10.
------------------	--


<b>SDLGV2188</b>	Green to Violet Loflosorb® 1kg Drum. Box of 10.
------------------	---

01491 527027

## Zero Dead Space Circuit for Rodents and Exotics

### Features & Benefits

- Delivers inlet gas directly into patient mask (zero dead space) for more rapid and predictable anaesthesia
- Reliable control of anaesthetic delivered to rodents
- Diaphragm openings available to suit a wide range of animals, from rats and mice to small animals and exotics
- More efficient single animal induction of rodents and other small animals

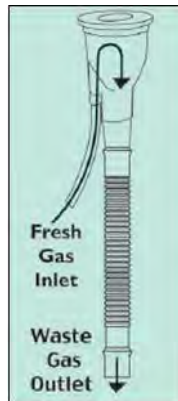

**Darvall ZDS (Zero Dead Space) Masks** provide predictable delivery of anaesthesia for rodents, birds and exotics 2kg and under. The stream-lined flow eliminates dead space, avoiding the need for tight fitting face masks.




Code	Description
ZDSM	Zero Dead Mask with 12mm and 20mm Diaphragms
ZDSDIA12	Replacement Diaphragm 12mm hole (to stretch to 15mm)
ZDSDIA20	Replacement Diaphragm 20mm hole (to stretch to 23mm)



## Exotic/Small Mammal Anaesthesia Mask and Non-Rebreathing Circuit

These especially small masks were specifically designed for use on small exotics such as pocket pets and birds. The small circuit hooks up the same way as a non-rebreathing circuit does to an anaesthesia machine. The masks fit tight around the nose and mouth. Soft conforming cone can be trimmed to custom fit animal's nose. The oxygen/gas tube position can be adjusted inside mask. Also circuit can be taped to the table to help keep it positioned on the animal's nose or muzzle.

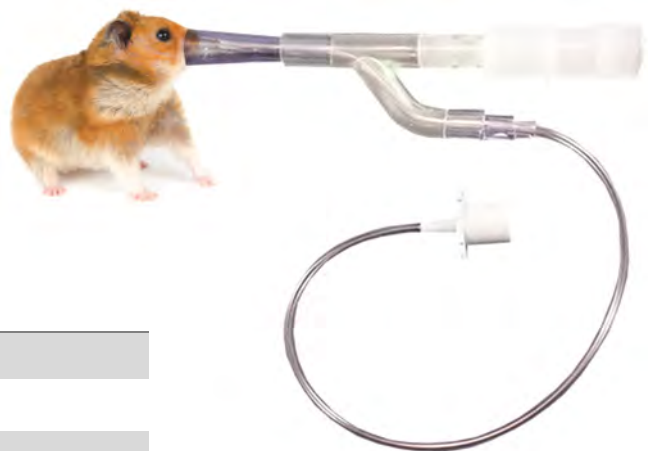

Code	Description
CIR-J0297AN	Exotic/ Small Mammal Circuit 9.0mm
CIR-J0297BN	Exotic/Small Mammal Circuit 12.5mm
CIR-J0297CN	Exotic/ Small Mammal Circuit 14.0mm

01491 527027

## In Circuit Pressure Manometer

- Works on common non-rebreathing circuits
- Simple to use – detach rebreathing bag and insert Manometer (directions included)

Code	Description
<b>PMD040</b>	In Circuit Pressure Manometer for 22mm Circuits. 0-60cm H2O
<b>PMD040TP</b>	In Circuit Pressure Manometer for Ayres T Piece. 0-60cm H2O
<b>PMD020</b>	In Circuit Pressure Manometer for 22mm Circuits. 0-25cm H2O
<b>PMD020TP</b>	In Circuit Pressure Manometer for Ayres T Piece. 0-25cm H2O

The In Circuit Patient Pressure Manometer is the safest way to accurately monitor the circuit pressure using non-rebreathers and circuits that do not utilize the manometer on the anaesthesia machine. It can also be used to add pressure monitoring capabilities to older anaesthesia machines that lack manometers.

The In Circuit Patient Pressure Manometers has a standard 22mm connector that can be used on all circuit types.

The Kit for the Ayres T Piece includes a connector to enable fitting to the smaller circuit.



## PEEP Valves

Intersurgical valves are 100% tested and individually identified with a code on the valve itself, so that performance data is available for future reference. The transparent body design, coupled with the coloured internal components, allow for easy blockage assessment and monitoring of patient's respiratory rate. Within the lightweight, robust L-shaped design we also feature an internal disc valve to satisfy the patient if the peak flow requirement exceeds the set flow.

Code	Description
<b>PEEP2202</b>	2.5cm H <sub>2</sub> O Colour Coded Black
<b>PEEP2205</b>	5.0cm H <sub>2</sub> O Colour Coded Blue
<b>PEEP2207</b>	7.5cm H <sub>2</sub> O Colour Coded Yellow
<b>PEEP2210</b>	10.0cm H <sub>2</sub> O Colour Coded Green
<b>PEEP2212</b>	12.5cm H <sub>2</sub> O Colour Coded Dark Pink
<b>PEEP2215</b>	15.0cm H <sub>2</sub> O Colour Coded Red
<b>PEEP2220</b>	20.0cm H <sub>2</sub> O Colour Coded Brown



01491 527027



## EMD Safety Valve



Most veterinary anaesthetic machines do not have a safety pressure relief valve. The pop-off valve can perform this function but it is easily defeated by being manually closed during mechanical ventilation and during system pressure testing. If the operator forgets to re-open the valve, tragedy results.

**The EMD Safety Valve is a simple solution to this problem.** It cannot be manually defeated. It is passive in that it does not need to be turned on like pressure alarms do. It can be placed anywhere in the rebreathing breathing circuit.

The valve opens once the airway pressure exceeds approximately 26 cm H<sub>2</sub>O, bleeding off gas pressure, thereby preventing barotrauma. Pressure will remain open until the pop-off valve is opened or the patient is disconnected from the circuit.

**All anaesthetic machines should be equipped with a Safety Valve.** Human error is unavoidable. This simple passive device will prevent a simple oversight from resulting in the death of a cherished pet and devastating effect a tragic medical error can have on the staff and the practice.

Code	Description
EMDSV	EMD Safety Valve
EMDSA	EMD Scavenging Adapter

## Intersurgical®

Single Code	Size	Box Qty	Box Code
INT2805	0.5L, 15mm Neck Breathing Bag for use with Ayres T Piece	10	INT2805/10
INT2810	1.0L, 22mm Neck Breathing Bag	10	INT2810/10
INT2811	1.0L, 15mm Neck Breathing Bag for use with Ayres T Piece	10	INT2811/10
INT2820	2.0L, 22mm Neck Breathing Bag	10	INT2820/10
INT2830	3.0L, 22mm Neck Breathing Bag	10	INT2830/10



Available individually or in box quantities.

## Latex Free Breathing Bag 0.5 litre with 22mm Neck

Code	Description
BBLF0522	0.5L with 22mm Neck



## Latex Free, Reusable Grey with Loop



Code	Size
BBLFR05	0.5L, 22mm Neck
BBLFR10	1.0L, 22mm Neck
BBLFR20	2.0L, 22mm Neck
BBLFR30	3.0L, 22mm Neck

Code	Size
BBBR05	0.5 litre with 22mm Neck
BBBR075	0.75 litre with 22mm Neck
BBBR1	1.0 litre with 22mm Neck
BBBR15	1.5 litre with 22mm Neck
BBBR2	2.0 litres with 22mm Neck
BBBR3	3.0 litres with 22mm Neck
BBBR4	4.0 litres with 22mm Neck
BBBR5	5.0 litres with 22mm Neck
BBBR6	6.0 litres with 22mm Neck
BBBR150	15.0 litres
BBBR30	30.0 litres

## Black Rubber



01491 527027

## Fibre Optic

### Right Hand Blades for Veterinary Use



Code	Size
<b>LBRH075</b>	75mm Blade
<b>LBRH100</b>	100mm Blade
<b>LBRH150</b>	150mm Blade
<b>LBRH200</b>	200mm Blade
<b>LBRH250</b>	250mm Blade
<b>LBRH300</b>	300mm Blade
<b>LBRH350</b>	350mm Blade
<b>LSETRH</b>	Set of 5 blades (75mm—250mm) and Standard Handle in a case

- Right Handed Veterinary Blades
- 75mm, 100mm, 150mm, 200mm, 250mm
- 300mm & 350mm blades also available
- Fibre Optic
- Choice of 3 Handles
- Autoclavable to 138°C

### Standard Fibre Optic (Human Type)

- 4000 Autoclave Cycle
- Compliant with Green System Handles

Code	Size
<b>LBFORMILO0</b>	Size 00
<b>LBFORMILO</b>	Size 0
<b>LBFORMIL1</b>	Size 1
<b>LBFORMIL2</b>	Size 2
<b>LBFORMIL3</b>	Size 3
<b>LBFORMIL4</b>	Size 4

Code	Size
<b>LHSL</b>	Slimline Handle
<b>LHSTDAU</b>	Standard Handle Autoclavable
<b>LHSTU</b>	Stubby Handle
<b>LHDPB</b>	BritePro Standard Handle



01491 527027

## BritePro Single Use Fibre Optic Blade & Handle Complete

Code	Size
LDPBMI00	BritePro Blade & Handle Miller 00
LDPBMI0	BritePro Blade & Handle Miller 0
LDPBMIL1	BritePro Blade & Handle Miller 1
LDPBMIL2	BritePro Blade & Handle Miller 2
LDPBMIL3	BritePro Blade & Handle Miller 3



Code	Size
LHDPB	BritePro Handle

BritePro Solo features a Single Use handle and LED light source with a fibre optic all-metal blade. The blade and the handle conform to all standard fibre optic systems.

## Conventional Laryngoscopes

### Miller Blades

Conventional Laryngoscopes have a physical bulb located in the blade.

Code	Size
LBCMIL0	Size 0, 75mm Blade
LBCMIL1	Size 1, 100mm Blade
LBCMIL2	Size 2, 155mm Blade
LBCMIL3	Size 3, 190mm Blade
LBCMIL4	Size 4, 200mm Blade



### Handles & Bulbs



Code	Size
LHCSL	Slimline Handle
LHCSTD	Standard Handle
BULBLG	Large Bulb
BULBSM	Small Bulb

01491 527027

## MRI Compatible Disposable Laryngoscopes

- Disposable system eliminates the risk of cross infection between patients
- Supplied with 2 x limited reuse 3.7V batteries, each providing 100 minutes of illumination
- Super bright LED far exceeds the light output requirements of ISO 7376
- ISO 7376:2020 Compatible, to be used with PROACT MRHPC blades only
- Tested to ASTM standards as suitable (MR Conditional) for use in MRI environments as a blade and handle set
- No touch removal and insertion of batteries reduces the risk of cross-contamination
- Non-magnetic, full metal handle construction with grooves offering strength and secure grip
- Affordable, cost effective, simple and disposable solution
- Supplied in easy to open, clearly identified MR sachets and dispenser box
- Latex Free

Available individually or in boxes of 10.

Single Code	Size	Box Qty	Box Code
<b>MRHPCMIL00</b>	Miller 00 Blade	10	<b>MRHPCMIL00/10</b>
<b>MRHPCMIL0</b>	Miller 0 Blade	10	<b>MRHPCMIL0/10</b>
<b>MRHPCMIL1</b>	Miller 1 Blade	10	<b>MRHPCMIL1/10</b>
<b>MRHPCMIL2</b>	Miller 2 Blade	10	<b>MRHPCMIL2/10</b>
<b>MRHPCMIL3</b>	Miller 3 Blade	10	<b>MRHPCMIL3/10</b>
<b>MRGSHB01</b>	MRI Compatible Handle Disposable with 2 x Battery	10	<b>MRGSHB01/10</b>
<b>MRBATT01</b>	Additional Battery Pack for MRI Handle		



01491 527027

## Portex Blue Line Ultra, Cuffed & Inner Cannula

These Portex Tracheostomy Tubes come complete with an Inner Cannula which are also available to purchase separately, in Packs of 2.

Code	Size
<b>TTCBLU060</b>	6.0mm
<b>TTCBLU070</b>	7.0mm
<b>TTCBLU075</b>	7.5mm
<b>TTCBLU080</b>	8.0mm
<b>TTCBLU085</b>	8.5mm
<b>TTCBLU090</b>	9.0mm
<b>TTCBLU100</b>	10.0mm



### Inner Cannula

Code	Size
<b>TTIC060</b>	6.0mm Pack of 2 Inner Cannula
<b>TTIC070</b>	7.0mm Pack of 2 Inner Cannula
<b>TTIC075</b>	7.5mm Pack of 2 Inner Cannula
<b>TTIC080</b>	8.0mm Pack of 2 Inner Cannula
<b>TTIC085</b>	8.5mm Pack of 2 Inner Cannula
<b>TTIC090</b>	9.0mm Pack of 2 Inner Cannula
<b>TTIC100</b>	10.0mm Pack of 2 Inner Cannula

## Jorvet, Silicone

- Most tubes are supplied with an obturator. These aid in insertion and are removed after placement.
- All tubes have a 15mm connector. This is the same fitting as an endotracheal tube and allows hook-up to an anaesthesia hose.
- All tubes have pliable wings or flanges allowing tying the tube in place with cotton (umbilical) tape or gauze.
- Tubes are non-cuffed or cuffed (needed when assisting ventilation).
- Tubes are single lumen cannula.
- All tubes are 100 percent silicone that softens at body temperature.
- All cuffed tubes are high volume low pressure and profile.

NOTE: Tracheostomy tubes are sized by their I.D. Tubes should be ~50% of the tracheal lumen width. The trachea changes diameter during the respiratory cycle. The tracheal body is wider than the larynx.

### Cuffed

Code	Size	Code	Size
<b>TTCS035</b>	3.5mm	<b>TTCS060</b>	6.0mm
<b>TTCS040</b>	4.0mm	<b>TTCS070</b>	7.0mm
<b>TTCS045</b>	4.5mm	<b>TTCS080</b>	8.0mm
<b>TTCS050</b>	5.0mm	<b>TTCS090</b>	9.0mm
<b>TTCS055</b>	5.5mm		



### Uncuffed

Code	Size	Code	Size
<b>TTUS025</b>	2.5mm	<b>TTUS055</b>	5.5mm
<b>TTUS030</b>	3.0mm	<b>TTUS060</b>	6.0mm
<b>TTUS035</b>	3.5mm	<b>TTUS070</b>	7.0mm
<b>TTUS040</b>	4.0mm	<b>TTUS080</b>	8.0mm
<b>TTUS045</b>	4.5mm	<b>TTUS090</b>	9.0mm
<b>TTUS050</b>	5.0mm	<b>TTUS100</b>	10.0mm

01491 527027

## Equine, Silicone Tracheostomy

Silicone softens at body temperature for patient comfort. Balloon inflates via one-way valve. The universal attachment fits large animal anaesthesia machine and ventilator hoses. Has a wide outer flexible plate for secure attachment

Code	Size
TTSE15	15mm Cuffed
TTSE20	20mm Cuffed
TTSE25	25mm Cuffed



## Metal Tracheostomy

Nickel-Plated Brass

Code	Size
TTNM18	18mm I.D. 22mm O.D.
TTNM22	22mm I.D. 28mm O.D.



Pape's Model, Self retaining.

## Tracheostomy Cleaning Brushes

Code	Size
TTCB080	8.0mm x 90mm
TTCB100	10.0mm x 90mm
TTCB120	12.0mm x 90mm
TTCB140	14.0mm x 90mm



## Hydrotrach HME Filter

The InterTherm T+ is a heat and moisture exchanger (HME) designed for use with spontaneously breathing patients in order to reduce loss of heat and moisture during respiration, ideal for patients and clinicians caring for patients living with a tracheostomy.

Small and lightweight, it features both a clipped suction port and swivel oxygen stem combined with a unique white corrugated paper HME design to provide the perfect combination between humidification output and low resistance to flow. The product is also available with supplemental oxygen tube.

Single Code	Box Qty	Box Code
INT1875	25	INT1875/25



01491 527027

## Surgivet Induction Chambers



**Code**      **Description**

**IC-V711801** Large Induction Chamber 4 x 4 x 12"

**IC-V711802** Small Induction Chamber 4 x 4 x 6"

## Economy

- Moulded one piece construction
- For induction, oxygen therapy or euthanasia
- Simple, secure closures
- Leak proof
- Easy to clean
- Internal Dimensions: 38cm x 24cm x 18cm.



**Code**      **Description**

**ICESM**      Economy Small



We can make Induction Boxes to the size of your choice. Please let us know your requirements by calling 01491 527027 or emailing on [sales@acevetssupplies.co.uk](mailto:sales@acevetssupplies.co.uk)



**01491 527027**



**INT1568** Straight connector 22F - 6mm oxygen stem

Available each or Box of 40



**INT1921** One-way directional valve 22M (flow in) - 22F (flow out)

Available each or Box of 75



**INT1945** Straight connector 15M-15M

Available each or Box of 30



**INT1948** Straight Connector 15M- 6mm swivel-15F

Available each or Box of 50



**INT1950** One-way directional valve 22M (flow in) - 22M (flow out)

Available each or Box of 40



**INT1961** Straight connector 22M-22F

Available each or Box of 40



**INT1963** Straight connector 22M - 6mm OD stem and cap - 22F

Available each or Box of 25



**INT1661** Straight connector 15M/11.5-10M. Fits into Ayres T Piece

Available each a Box of 35



**INT1943** Straight connector 15M - 22M

Available each or Box of 40



**INT1946** Straight connector 15M-15F

Available each or Box of 30



**INT1949** Straight connector 15F - 15F

£1.30 each

£26.00 box of 25



**INT1960** Straight connector 22M-22M

Available each or Box of 40



**INT1962** Straight connector 15M - 22F

Available each or Box of 30



01491 527027

## 40 Connectors– Straight

**INT1964** Straight connector 22M - 22F with 7.6mm Port and Cap

Available each or box of 25



**INT1966** Straight connector 22M - 7.6mm Port and Cap - 22M

Available each or box of 25



**INT1968** Tube Connector 22M with 6mm Oxygen Stem

Available each or Box of 40



**INT1970** Straight Scavenging Connector 22M-30M

Available each or Box of 35



**INT1971** Straight Scavenging Connector 22M-30F

Available each or Box of 25



**INT1977** Straight Connector 22M - 10mm swivel elbow - 22F

Available each or Box of 35



**GE2105488-016** Straight connector 22M - 22F with Port

Available each or Box of 10



**INT1965** Straight connector 22M - 6mm OD stem and cap - 22M

Available each or box of 25



**INT1967** Straight connector 22F-22F

Available each or Box of 35



**INT1969** Straight connector 22M/22M-15F

Available each or Box of 40



**INT1971001** Straight Scavenging Connector 30M-30F

Available each or Box of 75



**INT1974** Straight Connector 22M/15F with large Oxygen Stem

Available each or Box of 45



**INT2713** Straight connector 15M - luer lock port with retained cap - 22M/15F

Available each or Box of 25



01491 527027

**INT1997** Fixed elbow 15M-22M/15F

Available each or Box of 30



**INT1891** 22F-22F Elbow Connector with 7.6mm Port.

Available each or Box of 25



**INT1892** 22M-22M Elbow Connector with 7.6mm Port.

Available each or Box of 25



**INT1897** Fixed elbow 15M-22M/15F with Double Swivel

Available each or Box of 30



**INT2714** Fixed elbow 15M - luer port - 22M/15F

Available each or Box of 25



**INT1898** Double swivel elbow 15M - double flip top cap with seal - 22M/15F

Available each or Box of 30



## Luer Connectors

**INT2710** Luer Port Elbow 7.6mm Port.

Available each or Box of 25




**LUERFF** Female to Female Luer



**PFLA** Female Luer Lock to Tapered Catheter End Fitting 10-20FrID

**PFLAL** Female Luer Lock to Tapered Catheter End Fitting 20-40FrID

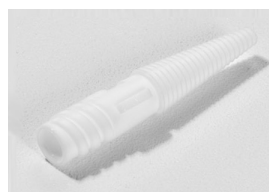

01491 527027

**INT1980** T-Piece 22M-22M-22M/15F

Available each  
or Box of 25



**INT1982** T-Piece 22M-22M-22F

Available each or Box of 20



**INT1902** 15mm paediatric traditional Y-piece 15F patient connection - 15M-15M

Available each or Box of 25



**YPP15M** 15mm paediatric traditional Y-piece 15F patient connection - 15M-15M

15M-15M

Available each or Box of 25



**INT1926** 22mm parallel Y-piece 22M/15F patient connection - 22M-22M

Available each or Box of 20



**INT1981** T-Piece 22M-22F-22M

Available each or Box of 25



**YPA22MLLP** 22mm anaesthetic Y-piece 22M/15F patient connection with luer lock port , Flexicare

Available each



**INT1903** 15mm paediatric traditional Y-piece 15F patient connection - 15M-15M with 7.6mm Ports

Available each or Box of 25



**INT1929** 22mm adjustable anaesthetic Y-piece with swivel neck 22M/15F patient connection - swivel 22M - swivel 22M with luer lock port

Available each or Box of 15



**OLC-Y** Oxygen Link Y connector 7mm—8mm OD. Pack of 10



01491 527027

## Needles for Blocks Under Ultrasound: 30° BEVEL (USB)

The USB stimulation needle features a familiar 30 degree 'back-cut' bevel tip for extra confidence in those used to this type of bevel, and improved perception of the needle's passage through tissues.

Code	Description	Diameter	Length
<b>USB050-22/10</b>	Ultrasound Needle 50mm. Box of 10	22G	50mm
<b>USB070-22/10</b>	Ultrasound Needle 70mm. Box of 10	22G	70mm
<b>USB090-21/10</b>	Ultrasound Needle 90mm. Box of 10	21G	90mm
<b>USB120-21/10</b>	Ultrasound Needle 120mm. Box of 10	21G	120mm

## Needles for Blocks Under Ultrasound: 20° BEVEL (USC)

The USC stimulation needle features a flattened 20 degree bevel tip, maximising the visibility under ultrasound – even when 'out of plane'. The needle is treated with a 'roughening' effect, maximising echographic properties and providing the user with a clear return.

Code	Description	Diameter	Length
<b>USC050-22/10</b>	Ultrasound Needle 50mm. Box of 10	22G	50mm
<b>USC070-22/10</b>	Ultrasound Needle 70mm. Box of 10	22G	70mm
<b>USC090-21/10</b>	Ultrasound Needle 90mm. Box of 10	21G	90mm
<b>USC120-21/10</b>	Ultrasound Needle 120mm. Box of 10	21G	120mm

## Needles for Blocks Under Ultrasound with Stimulation (Hybrid)

Code	Description	Diameter	Length
<b>EH035-22/10</b>	Hybrid Needle 35mm. Box of 10.	22G	35mm
<b>EH055-22/10</b>	Hybrid Needle 55mm. Box of 10.	22G	55mm
<b>EH095-21/10</b>	Hybrid Needle 95mm. Box of 10.	21G	95mm
<b>EH125-21/10</b>	Hybrid Needle 125mm. Box of 10.	21G	125mm
<b>EH150-21/10</b>	Hybrid Needle 150mm. Box of 10.	21G	150mm



01491 527027

## Echoplex Needles

Designed for clinicians who use Ultrasound and stimulation techniques to locate the peripheral nerve sheath.



Code	Size (G)	Length (mm)	Box Quantity
6194.103	21	100	20
6194.123	21	120	20
6194.153	21	150	20
6194.253	23	25	20
6194.353	23	35	20
6194.503	22	50	20
6194.853	21	85	20

## Echoplex Needles

Nerve stimulator allowing nerve location through the emission of a low intensity electrical current.



Code	Description
NS-750131	Plexygon Nerve Stimulator
08750302A	Replacement Cable (Red & Black)

## BSmart™ Injection Pressure Monitor

BSmart™ Injection Pressure Monitor is an easy-to-use inline injection pressure monitoring device providing clinicians with objective pressure information regardless of who performs the actual block.

While perineural injections are typically associated with low opening injection pressure (< 15psi) <sup>1</sup>, high opening pressure can indicate needle-nerve-contact <sup>1</sup> or intrafascicular needle placement.<sup>4</sup>

Monitoring the opening injection pressure with BSmart™ may prevent injections into poorly compliant tissue by aborting an injection and repositioning the needle when the monitor indicates that opening injection pressure is abnormally high (> 15psi).



Code	Description
BSMART	BSmart Injection Pressure Monitor

01491 527027

## Mistral Air Warming Unit

The Forced Air Warming unit is lightweight, easy to clean and made of shock absorbent material.

Improve patient outcome through safe, comfortable and smart patient warming.

Optimised heat transfer due to Multi Sensor Technology & Powerful Fan

Optimal elimination of air contaminants by the HEPA filter

Ergonomic and flexible in use as a result of the slim & lightweight design

Equipped with a versatile mounting bracket, the slim warming unit offers caregivers a diversity to securely position the warming unit across all care areas.




Code	Description
WAS001	Mistrial Warm Air Heater
WAS-TM	Trolley Mount for Mistrial
KK-J1304Q	Kennel Connector for Mistrial Warm Air System

01491 527027

## BAJA Disposable Blankets

PVS Order #78909274

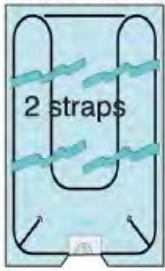


**Tube (M)**  
2-401 25" x 35"

PVS Order #78907456




**Tube (S)**  
2-400 20" x 25"

PVS Order #78907455

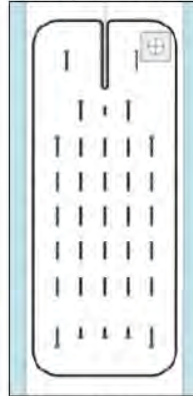

**Over/Under (L)**  
2-304 30" x 63"

PVS Order #78907454

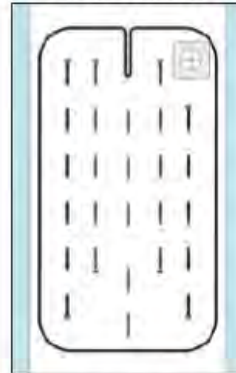

**Over/Under (M)**  
2-303 30" x 48"

PVS Order #78907453



**Over/Under (S)**  
2-302 30" x 31"



Code	Size
<b>WABBTM</b>	Medium Tube Under Blanket 635mm x 890mm
<b>WABBTS</b>	Small Tube Under Blanket 510mm x 635mm
<b>WABBOUL</b>	Large Over/Under Blanket 760mm x 1600mm
<b>WABBOUM</b>	Medium Over/Under Blanket 760mm x 1220mm
<b>WABBOUS</b>	Small Over/Under Blanket 760mm x 790mm

## BAJA Accessories Kennel Konnector



The BAJA Kennel Konnector is your conduit to safe patient warming inside the kennel.

- The blower and hose stay outside the Kennel
- The warming blanket, under/or around the patient is inside
- Using quick release reinforced rubber bands and foam inserts, the 2 piece door clamp assembly attached through the kennel door
- The 65mm hose (included) delivers air to the blanket, through the

Code	Size
<b>KKBPWS</b>	Kennel Konnector for the Baja Patient Warming System
<b>KKRS</b>	Kennel Konnector Replacement Straps
<b>KKRF</b>	Kennel Konnector Replacement Foam (Pair)
<b>KKDCA</b>	Kennel Konnector Replacement Door Clamp Assembly



Door Clamp Assembly



Foam

01491 527027



## Jorvet Reusable Warming Blankets

PART NO	SIZE
<b>WAB12060</b>	120cm x 60cm
<b>WAB9060</b>	90cm x 60cm
<b>WAB7570</b>	75cm x 70cm
<b>WABSET</b>	Set of 3 as above



The Jorvet blankets are made of a silky polyester fabric with an inlay of nylon water proof material to keep the air-flow circulating upward towards just the patient. These blankets are machine washable (hang to dry) therefore much more economical than other air blankets on the market, which are disposable. There is a 5" diameter round opening with a nylon cord to secure the blanket to any of the major brands of warm air convection units.

## Hoverheat Forced Air Warming

HoverHeat is the only forced air warming system that provides true under-body warming. With its proprietary internal components, HoverHeat provides a cushion of heated air to warm under and over body - something not previously possible with forced-air warming!

The most Complete patient air warming system ever  
True under and over body warming with one warm air blower

Use any warm air blower

Use for surgical, dental, endoscopy, and imaging procedures

Use for Pre-op, Intra-op, and Post-op warming



### The set includes:

- 1 Small HoverHeat.
- 1 Medium HoverHeat.
- 1 Large HoverHeat.
- 1 HoverHeat Connector.

Code	Description
<b>WABHH2914</b>	Hoverheat Warming 29" x 14" Small
<b>WABHH3517</b>	Hoverheat Warming 35" x 17" Medium
<b>WABHH4020</b>	Hoverheat Warming 40" x 20" Large
<b>HHCON</b>	Hoverheat Connector
<b>HHADA</b>	Hoverheat Adapter
<b>HHSET</b>	Hoverheat Set



01491 527027

## Foil Blankets

PART NO	DESCRIPTION
FBM	130cm x 210cm
FBS	120cm x 90cm
WBHP-NB	Wheat Bag, Rectangle –Navy Blue



## Reflect-a-Therm Table Cover

- Reflect-a-Therm Table cover offers excellent insulation for the animal.
- The reflective layer and the insulating bubbles stops up to 80% of the animals heat escaping through the cold stainless steel.
- Excellent way of keeping animals warm during examination, where Hyperthermia is a real risk.
- The PVC is sealed using HF welding to ensure the cover is claw proof and leak proof.
- Can be easily wiped clean and is padded for added comfort.



PART NO	SIZE
RATTC11560	115cm x 60cm
RATTC12050	120cm x 50cm

## Silver Heat Products

Looking for an **inexpensive, easy-to-use method of improving patient warming**, in any situation, **Silver Heat products** are made from a special, aluminium-coated reflective nonwoven material that bounces up to 80% of your patient's radiant heat loss back towards the animal, providing what we like to call "patient-powered" warming.

### Leggings



Code	Size
SHLXL	Legging, Extra Large 760mm. Pack of 2
SHLL	Legging, Large 560mm. Pack of 2
SHLM	Legging, Medium 380mm. Pack of 2
SHLS	Legging, Small 280mm. Pack of 2
SHLXS	Legging, Extra Small. 180mm. Pack of 2
SHPAP	Papoose, 330mm x 610mm
SHNEST	Small Animal Nest, 310mm dia

### Papoose



- Excellent for patients too small to cover with a blanket
- Velcro closures
- Adjustable sizing yields good fit for patients under 10kg
- Great for shorter dental procedures

### Small Animal Nest



- Warm, Cosy Comfort for Small Cats and Exotics
- Circular configuration accommodates small animals of all shapes
- Elastic closure at top ensures closure around animal
- Full body radiant heat reflection

01491 527027

## Thermal Mattress

For Surgery, Endoscopy, Dental, Imaging  
Radiolucent X-Ray and CT  
Adjunct to Active Warming  
Easily Cleaned  
Reusable

- Available in 3 Sizes
- Small 510mm x 710mm
- Medium 610mm x 860mm
- Large 710mm x 1015mm

Two Layer Internal Design  
Insulative conductive thermal barrier  
Radiant Heat Protective Barrier



Code	Size
<b>WTB-LAR</b>	Warm Thermal Blanket Large 710mm x 1015mm
<b>WTB-MED</b>	Warm Thermal Blanket Medium 610mm x 860mm
<b>WTB-SML</b>	Warm Thermal Blanket Small 510mm x 710mm
<b>WTB-SET</b>	Set of Small, Medium & Large Blankets

- Use either ConRad Thermal Blanket immediately upon induction of anaesthesia or sedation.
- Prevents the rapid Phase I drop in body temperature.
- The time taken to prep for surgery and apply monitors can drop the animal's body temperature several degrees even before surgery or imaging starts, unless preventative measures are taken.
- The General Use ConRad Thermal Blanket is recommended for all procedures and imaging, except for MRI.
- The MRI-Safe ConRad Thermal Blanket (Page 62) is used with MRI procedures.

## MRI Safe Thermal Mattress

Designed exclusively for use in MRI  
No MRI Signal Loss  
Clinically Tested to Maintain Body Temperature  
Available in 3 Sizes

- Small 510mm x 710mm
- Medium 610mm x 860mm
- Large 710mm x 1015mm

Two Layer Internal Design  
Insulative conductive thermal barrier  
Proprietary MRI Safe heat loss barrier



Code	Size
<b>WABMRI-LAR</b>	Warm Thermal Mattress Large 710mm x 1015mm
<b>WABMRI-MED</b>	Warm Thermal Mattress Medium 610mm x 860mm
<b>WABMRI-SML</b>	Warm Thermal Mattress Small 510mm x 710mm
<b>WABMRI-SET</b>	Set of Small, Medium & Large Mattresses

01491 527027

## Dormosafe Electric Blanket

The Dormosafe heating mat is the latest development of professional heating systems for veterinary medicine. The padded mat combines safety, user-friendliness and heating quality. Three temperature sensors distributed throughout the mat continuously measure the accumulated heat between the patient and the mat. The redundantly protected microcontroller in the mat regulates the heat output depending on the degree of cooling of the patient. The pre-set temperature of 42°C is thus neither exceeded nor fallen short of, which guarantees an even and safe heating.

- Safety Voltage—electrical safety for patient and personal.
- Overheat protection with 3 Temperature Sensors—optimise heat transfer to 42°C.
- Acoustice Alarm—due to errors or temperature above 45°C.
- Size 45cm x 35cm x 2cm



Code	Size
DHM45352	Dormosafe Heat Mat Complete 45cm x 35cm x 2cm
DHM60452	Dormosafe Heat Mat Complete 60cm x 45cm x 2cm

### Warm Patients Heal Faster

A constant body core temperature is a basic prerequisite for the function of vital organs. It reduces the risk of surgical complications, accelerates the recovery phase and healing.

### Easy, Fail-Safe Handling

The control electronics of the Dormosafe heating mat are accommodated in the mat itself so there is no longer any need for an external control unit. The flexible mat is therefore completely insensitive to falls. Due to the three sensors, a wrong positioning of the patient is impossible and therefore handling is easier.

## Heated Veterinary Operating Tables

A simple but effective product that gently warms small mammals on the operating table minimising loss of body heat, and increasing survival rates.

The V500DVstat is constructed from tough ABS plastic with an anodised aluminium panel bonded to the surface to aid cleaning and sterilisation. A separate digital control module allows the precise temperature of the surface to be selected by the surgeon within the range 10 to 50 degrees C, and operation of the heater is indicated by an amber light.

Size: 350mm x 530mm

The V500 is constructed from tough ABS plastic, with an anodised aluminium panel bonded to the surface to aid cleaning and sterilisation. This is heated and controlled by an electronic thermostat to 35 degrees C. An amber light displays operation of the heater.




Code	Description
V500VSTAT	Heated Table 350 x 530mm with Digital Control
V500	Heated Table 350 x 530mm to 35°C

01491 527027

## HUG-U-VAC

- Conforms Gently to Body Contours
- Pressure Sore Prevention
- Body Heat Conservation
- Emergency Fracture Stabilization
- Secure Head and Neck Support
- X-Ray Lucent
- Easily Cleaned
- Pre-op Prep – shaving or sterile scrub
- Surgery – ventral, dorsal or lateral position
- Post-op Recovery – prevent aspiration
- Radiological Exams – precise positioning for accurate x-rays
- Ultrasound and CT scans
- Emergency Transport
- Newborn Care – positioning and heat conservation

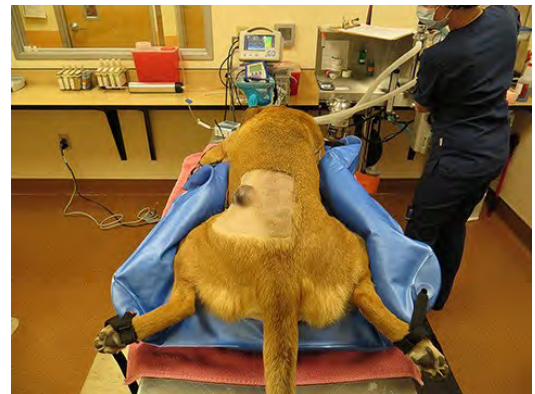

Code	Description
<b>HUV1624</b>	Small Positioner 16" x 24" Puppies & Kittens
<b>HUV2030</b>	Medium Positioner 20" x 30" Cats & Small Dogs
<b>HUV2436</b>	Large Positioner 24" x 36" Medium and Large Dogs, Goats, Sheep
<b>HUV3045</b>	Extra Large Positioner 30" x 45" Large Dogs, Goats & Sheep
<b>HUV3654</b>	Jumbo Positioner 36" x 54" Large Animals, Alpacas, Llamas, Fools & Mini Horses
<b>HUVHHP</b>	Head Positioner
<b>HUVMVP</b>	Manual Vacuum Pump & Hose
<b>HUVVDB</b>	Duffel Bag
<b>HUVPK</b>	Patch Kit
<b>HUVELRS</b>	Extra Limb Restraints (Set of 4)
<b>HUVSET</b>	5 Body Positioners, 1 Head Positioner, 1 Vacuum Pump, 1 Duffel Bag, 1 Patch Kit and 20 Limb Restraints



01491 527027

## Buster Vacu Support

Ideal support for surgery and examination of patient. Keeps the patient insulated from cold operating tables. Resistant to all commonly used disinfectants. Radiolucent.

Code	Description
BB3060	30cm x 60cm
BB50100	50cm x 100cm
BBVP	Vacuum Pump (Yellow)
BBVP1	Vacuum Pump (Blue)



Foot Pump (Blue)



Foot Pump (Yellow)

## Schure Slide™ Patient Transfer System

Used to transfer patients safely and easily. The sliding action eliminates the need to lift the patient, removing the risk of possible back injury to personnel. There are no steel rollers or hard surfaces to cause pain or discomfort to the patient, and the unique patented design allows patient transfer in any position. The Schure Slide™ is radiolucent and is non-staining.

**Measures:** 46 x 84 cm

Code	Size
PTS4684	Schure Slide Patient Transfer System



## Redmark Cradles

Small animal cradle system for X-Ray and Surgical Positioning.

Code	Description
CR5050	Large Animals over 32kg weight. 50cm x 50cm
CR4540	Medium Animals 13kg—32kg weight. 45cm x 40cm
CR3427	Small Animals up to 13kg weight. 34cm x 27cm
CRSET3	Set of 3 Small, Medium, Large Cradles



01491 527027

## Oxygen Link Tubing

Code	Description
<b>OLT30M</b>	Oxygen Link Tubing 30m
<b>OT1180</b>	Oxygen thick wall tube, 5mm ID - 8mm OD



OLT30M

The 3.0mm inner diameter tubing comes in a dispensing box with a 30m coil that can be cut to the desired length.



OT1180

1 roll of 50m length.

## Nasal Cannulas

Code	Size
<b>NC-LGE</b>	Nasal Cannula—Large (Adult)
<b>NC-MED</b>	Nasal Cannula—Medium (Paediatric)
<b>NC-SML</b>	Nasal Cannula—Small (Infant)
<b>NC-XSM</b>	Nasal- Cannula—Extra Small (Neonatal)



## Oxygen Link Tubing Connector

Code	Description
<b>OLC-Y</b>	Oxygen Link Tubing. Connector. Pack of 10.



## Nasal Oxygen Catheter

Code	Size
<b>NOC0639</b>	6fr x 39cm
<b>NOC0839</b>	8fr x 39cm
<b>NOC1039</b>	10fr x 39cm



01491 527027

## Bubble Humidifier Bottle

Code	Description
<b>OHR0788</b>	Oxygen Humidifier Bottle Disposable
<b>OHB21135</b>	Oxygen Humidifier Bottle Autoclavable



Autoclavable

Disposable

## Oxygen Flowmeters

Code	Description
<b>OFDP015</b>	Oxygen Flowmeter 0 - 15LPM. Direct Probe with Tube Nipple
<b>OFBSP015</b>	Oxygen Flowmeter 0 - 15LPM. 3/8" BSP with Tube Nipple



Direct Probe

3/8" BSP

## Regulators

Code	Description
<b>O2RBNBSP</b>	Oxygen Regulator Bullnose with 3/8" BSP



## Tubing Nipple

Code	Description
<b>O2TN</b>	Tubing Nipple



## O2 Splitters

Code	Description
<b>O2TWIN</b>	Twin Schrader Adapter + BS Probe
<b>O2TWIN-MS</b>	Schrader/Mini Schrader + BS Probe
<b>O2TWIN-MS2</b>	Mini Schrader Adapter + BS Probe



O2TWIN



O2TWIN-MS



O2TWIN-MS2

01491 527027



## Buster ICU Cage

Buster ICU Cages are collapsible, portable and can be set up in the prep room or any appropriate sit.

They fit standard stainless steel kennels.

They allow for the provision of oxygen or heat and also allow for the cooling of animals.

They can also be used as an anaesthetic chamber provided a co-axial circuit is used.

### Includes...

1 x Carry Bag

1 x ICU Cage

1 x Mat

Buster Gel Packs (1 x S, 2 x M, 4 x L)

Oxygen Accessories Pack

6 x Venturi Nebuliser, T Piece, Cap 2.1m oxygen tubing.

BUSTER Electric heat pad—Medium and Large BUSTER ICU Cages only.

See Page 89 for the Vetronic Oxygen Analyser

Code	Size
<b>ICU-LGE</b>	ICU Cage Large
<b>ICU-LGEHM</b>	ICU Cage Large with Heat Mat
<b>ICU-MED</b>	ICU Cage Medium
<b>ICU-MEDHM</b>	ICU Cage Medium with Heat Mat
<b>ICU-SML</b>	ICU Cage Small



## Bodok Seals & 'O' Rings

Code	Description
<b>BODOK</b>	Bodok Seal
<b>ORINGBN</b>	Bullnose O Ring
<b>ORINGVAP</b>	Vaporiser O Ring



## Cylinder Keys

Code	Description
<b>CKBN</b>	For Bullnose Cylinders
<b>CKPI</b>	For Pin Index Cylinder (Ratchet Type)
<b>CKPIE</b>	For Pin Index Cylinder (Economy)



Bullnose



Pin Index Ratchet



Pin index Economy

## Jugular Catheters—Vygon

### Description

Multicath Expert is the expert silver ion impregnated catheter. It possesses benefits such as:

- Pressure-rated
- Long-lasting; effective up to 30 days
- Atraumatic tip meaning its soft, flexible design allows smooth transition over the guidewire whilst reducing the potential for internal vessel trauma
- Anti-kink Teflon coated guidewire

Multistar antibiotic technology; releases rifampicin and miconazole on both inner and outer surfaces at a slow, steady rate.

Code	Lumen	Length (cm)	Size (Fr)	OD (mm)
<b>CVC1202.063</b>	2	6	4.5	1.5
<b>CVC1203.082</b>	2	8	4.5	1.5
<b>CVC1203.123</b>	2	12.5	4.5	1.5
<b>CVC157.508</b>	2	8	5.5	1.85
<b>CVC1233.061</b>	3	6	4.5	1.5
<b>CVC1233.081</b>	3	8	4.5	1.5
<b>CVC1233.101</b>	3	10	4.5	1.5
<b>CVC1231.121</b>	3	12.5	4.5	1.5
<b>CVC155.167</b>	3	16	7.5	2.5
<b>CVC155.207</b>	3	20	7.5	2.7
<b>CVC0155.508</b>	3	8	5.5	1.8
<b>CVC0155.512</b>	3	12.5	5.5	1.8



### Contents::

- 1 x Multicath forward catheter
- 1 x Introducer needle 18G x 70mm
- 1 x Mosquito introducer cannula 18G x 64mm
- 1 x 'J' guidewire in feeder; length 60cm, OD 0.9mm
- 1 x Vessel dilator
- 2 x Latex-free injection membranes
- 1 x Mini scalpel N
- 1 x 5mL luer-slip syringe
- 1x secondary fixation wing and anchor

01491 527027

# Jugular Catheters / CVCs—Infusion Concepts

Kit contains: Catheter, movable suture wing, scalpel and handle, guide wire, Seldinger introducer, dilator, and non latex injection caps.

Code	Lumen	Size
JC1308	1	3FR (20G), 8cm
JC1310	1	3FR (20G), 10cm
JC1313	1	3FR (20G), 13cm
JC1320	1	3FR (20G), 20cm
JC1408	1	4FR (18G), 8cm
JC1410	1	4FR (18G), 10cm
JC1413	1	4FR (18G), 13cm
JC1515	1	5FR (16G), 15cm
JC1530	1	5FR (16G), 30cm
JC1620	1	6FR (14G), 20cm
JC2406	2	4FR (22&22G), 6cm
JC2408	2	4FR (22&22G), 8cm
JC2410	2	4FR (22&22G), 10cm
JC2412	2	4FR (22&22G), 12cm
JC2510	2	5FR (18&22G), 10cm
JC2515	2	5FR (18&22G), 15cm
JC2720	2	7FR (14&22G), 20cm
JC2730	2	7FR (14&29G), 30cm

Code	Lumen	Size
JC3510	3	5.5FR (20, 22 & 22G), 10cm
JC3515	3	5.5FR (20, 22 & 22G), 15cm
JC3715	3	7FR (16, 18 & 18G), 15cm
JC3720	3	7FR (16, 18 & 18G), 20cm
JC3730	3	7FR (16, 18 & 18G), 30cm



## Guide Wires

Code	Make	Size
GW350	Infusion Concepts	3F x 50cm (0.018")
GW460	Infusion Concepts	4F x 60cm (0.021")
GW550	Infusion Concepts	5F x 50cm (0.025")
GW6770	Infusion Concepts	6-7F x 70cm (0.035")
WG1170	Surgivet	0.017" x 70cm for use with 5fr or 6fr
WG3860	Surgivet	0.038" x 60cm for use with 8fr
WG4560	Surgivet	0.045" x 60cm for use with 10fr or 12fr
GW0923.50	Vygon	0.46mm x 50cm
GW0923.60	Vygon	0.9mm x 60cm



## Arterial Catheters

Code	Size
AC2080	Catheter set 20G, 80mm
AC2280	Catheter set 22G, 80mm



01491 527027

## Surgivet Clear View Foley Catheters

Closed-end silicone catheter with distal balloon cuff used for intermittent drainage of Urine from the bladder or contrast studies of urogenital tract.

- Ribbed balloon provides even inflation and aids in retaining the catheter inside the bladder.
- Recommended to only be used with PTFE coated wire guides
- Sterile
- Radiopaque


Available individually or in box quantities



Code	Size	Box Qty	Box Code
<b>FCS0530</b>	5F x 30cm	10	<b>FCS0530/10</b>
<b>FCS0555</b>	5F x 55cm	10	<b>FCS0555/10</b>
<b>FCS0630</b>	6F x 30cm	10	<b>FCS0630/10</b>
<b>FCS0655</b>	6F x 55cm	10	<b>FCS0655/10</b>
<b>FCS0830</b>	8F x 30cm	10	<b>FCS0830/10</b>
<b>FCS0855</b>	8F x 55cm	10	<b>FCS0855/10</b>
<b>FCS1030</b>	10F x 30cm	10	<b>FCS1030/10</b>
<b>FCS1055</b>	10F x 55cm	10	<b>FCS1055/10</b>
<b>FCS1230</b>	12F x 30cm	10	<b>FCS1230/10</b>
<b>FCS1255</b>	12F x 55cm	10	<b>FCS1255/10</b>
<b>FCS1455</b>	14F x 55cm	10	<b>FCS1455/10</b>
<b>FCS1655</b>	16F x 55cm	10	<b>FCS1655/10</b>
<b>FCS1855</b>	18F x 55cm	10	<b>FCS1855/10</b>
<b>FCS2065</b>	20F x 65cm	10	<b>FCS2065/10</b>
<b>FCS2265</b>	22F x 65cm	10	<b>FCS2265/10</b>

## Clear Silicone with Guide Wires

Code	Size
<b>FCS0630GW</b>	6F x 30cm
<b>FCS0655GW</b>	6F x 55cm
<b>FCS0830GW</b>	8F x 30cm
<b>FCS0855GW</b>	8F x 55cm
<b>FCS1030GW</b>	10F x 30cm
<b>FCS1055GW</b>	10F x 55cm
<b>FCS1235GW</b>	12F x 35cm

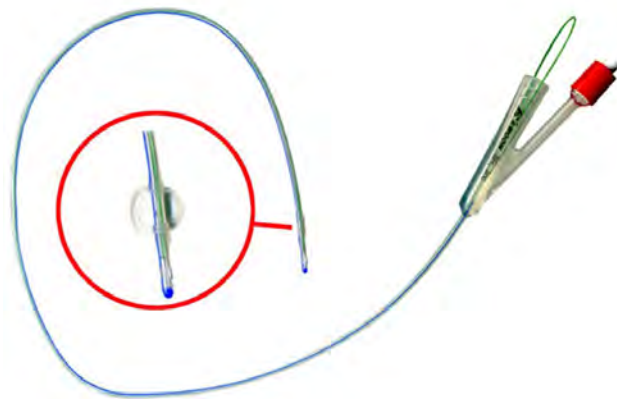

01491 527027

## Buster Dog Catheters

Female luer lock fitting

- Rigid though flexible polypropylene catheter tubing for easy introduction
- Rounded end and atraumatic side holes to avoid tissue damage
- User friendly blister packaging
- 12/pk



Code	Description
<b>CDU13500</b>	1.3mm x 500mm (4FG)
<b>CDU20500</b>	2.0mm x 500mm (6FG)
<b>CDU26500</b>	2.6mm x 500mm (8FG)
<b>CDU33500</b>	3.3mm x 500mm (10FG)

## KatKath™ Catheter

Code	Description
<b>KAT3511</b>	Standard KatKath 3.5FR x 11cm. BOX OF 10
<b>KAT3514</b>	Standard KatKath 3.5FR x 14cm BOX OF 10
<b>KAT3518</b>	Adjustable KatKath 3.5FR x 18cm. Box of 10



Standard

Simple to use—unblock the cat, leave in place

Non-reactive material—low friction. Vital when dealing with an inflamed oedematous urethra

End hole—allowing concurrent flushing during placement

Available with a rotating suture collar—allowing natural rotation when in situ

Kink and crush resistant without the use of a stylet

Tapered distal tip—allowing atraumatic catheter placement

Adjustable has an adjustable wing



Adjustable

## 'Little Herbert'

Allows for firm connection of drainage lines to urinary catheters without interfering with the comfort factor.

Code	Size
<b>TCLLA</b>	Little Herbert Luer Lock Adapter



## IV Pressure Bags

### Riester Matpak

- Sturdy cotton cuff in blue, with suspension loop, washable up to 60°C
- Precise reading of the level in the bag
- Chrome-plated manometer with easy-to-read 49 mm diameter aluminium scale up to 300mm Hg
- Specially hardened copper-beryllium membrane
- No zero point setting required
- Latex pad, resistant to pressure
- Latex-free inflation bulb



Code	Description
IVPB05	Riester Matpak 500ml
IVPB10	Riester Matpak 1000ml
IVPB30	Riester Matpak 3000ml
IVPB50	Riester Matpak 5000ml

### Semi Disposable

Code	Description
IVPBD05	Semi-Disposable 500ml
IVPBD10	Semi-Disposable 1000ml
IVPBD30	Semi-Disposable 3000ml

- Clear sleeve around bag enables users to view pressurised fluids and text
- Durable construction, will not rip or tear
- Palm size, high volume bulb for quicker inflation
- Sturdy loop hanger
- Latex free, Smooth cuff material offers wipe-clean surface
- Easy to read, fully enclosed pressure gauge enables 360 degree viewing



### Clear Re-usable

Code	Description
IVPB05CC	500ml Clear Pressure Bag
IVPB10CC	1000ml Clear Pressure Bag



### Mantis IV Hanger

Code	Description
MIVH	Mantis IV Hanger



### Extension Sets

Code	Description
ESM25200	Extension Set, Mini Bore 200cm
ESWB4125	Extension Set, Wide Bore 4.0mm x 1255mm



01491 527027

Luer-lock T connector with Bionector allows intermittent injection via latex-free membrane during simultaneous infusion. T connectors are ideal for general veterinary ward and theatre use as they can be pre-placed and heparinised by nursing staff on admission, even before a plan has been made regarding drug and infusion administration.

They are a universal IV solution, veterinary surgeons are able to administer drugs and rapidly connect infusion lines at will to any patient. Significant time savings are possible in veterinary wards and theatres if in-patients routinely have a Bionector T piece in place.

The Bionector needle-free device allows easy IV disconnection for exercise breaks, and can be disinfected with most antiseptic solutions (example: 70° alcohol). Avoid using antiseptic containing both alcohol and hydrogen peroxide.



## Single Bionector

Code	Description
0896.03	Bionector Closed System Needle Free Device



## Bionector T-Piece

Code	Description
1VTUSBN-448	10cm Luer-Lock T-Piece Extension Connector with Bionector
1VTUSBN-450	Bionector T-Piece with Two Needle-free Access Ports



## Bionector and Spike Bag

Code	Description
0856.11	Bionector and Spike Bag



Code 856.01: spike with Bionector, for aspirating or adding fluids to bags or bottles of intravenous solutions. Can also be used for paediatric bags.

The needless connector Bionector maintains at all times a closed system allowing to aspirate or adding fluids in a aseptic manner.

Scan here for more product info / video and brochure →



01491 527027

## Bionector Vial Caps

For use when aspirating drugs from capped vial and drug reconstitution.

Code	Size
0821.01	14mm for use with 5ml bottles
0821.02	20mm for use with 10—100ml bottles



## 3M™ CuroS™ Disinfecting Port

Disinfects in 1 minute. Protects I.V. access points for up to 7 days if not removed.

Helps reduce the risk of contaminants from entering the I.V. port through the catheter, post-insertion

Brightly coloured to help clinicians verify and measure compliance



Code	Description
WWJIM001-270	3M™ CuroS Disinfecting Port

## BD 3 Way Connecta™

Code	Description
WWJIM001-270	3m™ CuroS Disinfecting Port

Available singularly or in boxes of 100.



## V Kit 20G Sedation Pack

Code	Description
VKIT20	V Kit 20G Sedation Park

Includes Cannula, IV, Jelco 20g, Pink Dressing, Adhesive IV, Tegaderm, 5cm x 5.7cm  
 Injection Cap, Female Luer  
 Syringe, 5ml Luer Slip, Prefilled, Posi-Flush



01491 527027



## Biopsy Needles

### Carefusion Jamshidi Bone Marrow/Aspiration Needle



Code	Size
NB11	11G x 10cm
NB13	13G x 8.9cm
NB15	15G x 7.9cm
NB18	18G x 7.9mm

### Jorvet

These bone marrow needles allow for an actual biopsy versus the aspiration other needles provide. A special reamer is included to allow retrieval of actual tissue sample.



Code	Size
NBA08	08G x 4"
NBA11	11G x 4"
NBA13	13G X 4"

### Surgivet Vet-Core Biopsy Needle



Used for soft tissue core biopsy.

Code	Size
VC1409	14GA x 9cm
VC1415	14GA x 15cm
VC1609	16GA x 9cm
VC1615	16GA x 15cm
VC1809	18GA x 9cm

01491 527027

## PalmSat VET Pulse Oximeter

The 2500A VET handheld pulse oximeter is a universal tool with unmatched versatility and cost-effectiveness for all veterinary SpO<sub>2</sub> monitoring and research applications. Compact and easy to use, the 2500A VET is proven accurate for pulse rates up to 450 beats per minute (BPM) and is ideal for monitoring during surgeries and dental procedures. As the industry benchmark for durability and performance, Nonin veterinary pulse oximeters have gained the confidence of veterinary professionals worldwide.

Code	Description
<b>MON2500AV</b>	PalmSat 2500AV Pulse Oximeter with Lingual Sensor and Black Carry Case
<b>MON2500AV-02</b>	PalmSat 2500AV Pulse Oximeter with Transflectance Sensor and Black Carry Case

## Spare Sensors

Code	Description
<b>SP02LS-NON1M</b>	Nonin Lingual Sensor 1M Cable
<b>SP02TP-NON1M</b>	Nonin Transflectance Sensor 1M Cable
<b>SP02WS-NON1M</b>	Nonin Wrap Sensor 1M Cable
<b>2500AVCASE</b>	Black Carry Case for PalmSat
<b>2500C</b>	Charger Stand Kit



See also Nonin SpO<sub>2</sub> sensors on Page 68.



Charger Kit & Battery Pack



01491 527027

# Edan H100B Veterinary Pulse Oximeter

Standard System:

- Main Unit
  - LCD Screen
  - RTC display
  - SpO2 and Pulse Rate Measurement
  - Display Numeric and Waveform of SpO2 Simultaneously
  - Convenient AA Size for up to 48 Hours Working
  - Up to 100 IDs and 300 hours data storage capability for continuous monitoring
  - The latest 10 minutes trend graph and table of SpO2 and Pulse Rate can be reviewed in the screen
  - Compatible with Nellcor VETSAT Veterinary sensor
  - Pitch tone
- Standard Accessories:

- EDAN Veterinary Pulse Oximetry
- Sensor
- AA Size Battery
- Carrying Case
- User Manual
- VE-H100B Quick Reference Card



Code	Description
<b>MON-H100B</b>	H100B Pulse Oximeter
<b>CASE-H100B</b>	Carry Case for H100B
<b>CHAR-H100B</b>	Charger Stand for H100B
<b>SILI-H100B</b>	Silicone Boot for H100B
<b>SP02EC-01-13.21001</b>	Sp02 Extension Cable
<b>SP02EP-02.01.110498</b>	Sp02 Vet Sp02 Sensor



Carry Case



Charger Stand



Silicone Boot

## SpO2 Sensor Clips

Code	Description
<b>SP02CLIPS</b>	Oxygen Sensor Clips Set of 2, Large & Small
<b>MAOSCLIP-NC</b>	Oxygen Sensor Clip for Nellcor
<b>V1709</b>	'C' Clip replacement for use with Surgivet V1707



SP02CLIPS



MAOSCLIP-NC



V1709

## SpO2 Sensors—Surgivet

Compatible with Surgivet Monitors V3304, V3402, V3404, V6004, V9004, V9200, V9230, V8400D, V3395 and V1030.



**SP02S-V1700** This reflectance sensor can be applied to the ventral tail base of a patient and secured with a cohesive bandage.

**SP02S-V1703** The Universal 'Y' sensor is designed to be used on areas of thicker tissue. It is the standard clip used in most small animal practices.



**SP02S-V1707** The 'C' sensor is perfect for cats, dogs and horses. The sensor should be placed on the thigh, metatarsal, or metacarpal for smaller patients. The sensor should be placed over the Achilles tendon, tongue, prepuce, vulva, or through the toe webbing on large patients.

## Mindray

### SP02EC-MR

Interface Cable for Mindray and Passport Monitors. Equivalent to Mindray 010-20-42594. To be used with SP02EPSQ-MR




### SP02EC-IMEC

2.4 metre interface Cable for Mindray IMEC. To be used with SP02EPSQ-MR.




### SP02EPSQ-MR

Patient End Lingual Probe equivalent to Mindray 115-036222-00. Supplied with Small and Large Clip.

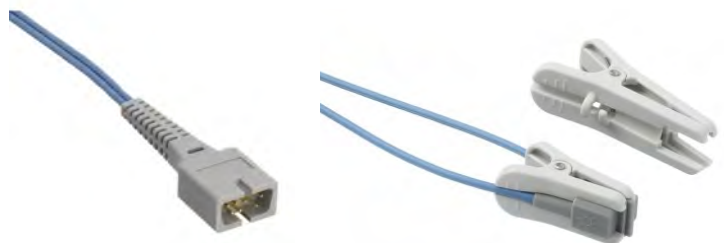

### SP02DC-MR

Direct Cable for Mindray and Passport Monitors. Supplied with Small and Large Clip.

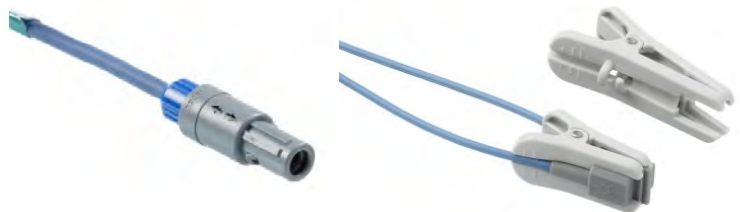

01491 527027

## Nellcor

### SP02-S810-01P0

Nellcor Direct Cable with Lingual Sensor



### SP02-S810-01P0

Nellcor Direct Cable with Lingual Sensor




### SP02EP-FD

Nellcor Lingual Sensor (Patient End)



## Veterinary Multi Site SpO2 Sensor

**SP02-10181 Nellcor** Multi Site Sensor Supplied with Small and Large CLIP.

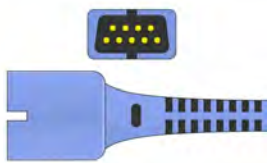

## Nonin

**SP02DC-NON09** Nonin Lingual Ear Probe. Supplied with Small and Large Clips. 90cm in length.



**SP02DC-NON30** Nonin Lingual Ear Probe. Supplied with Small and Large Clips. 3 metres in length.

## GE Healthcare/Datex

**SP02EPDC-DX** Direct 3 metre Cable equivalent to GE TS-E4-N. Supplied with Small and Large Clip



**SP02EPDC-GE** Direct 3 metre Cable equivalent to GE TS-E4-GE. Supplied with Small and Large Clip



**SP02TS-E4-N** Direct 4 metre Cable with Lingual Probe. *Genuine GE Part.*



**SP02TS-E4-GE** Direct 4 metre Cable with Lingual Probe. *Genuine GE Part.*



**SP02IC-NCDX** Interface Cable for GE Monitors with Nellcor Technology (use with SP02EP-FD on Page 68)



01491 527027

## Masimo



### SP02EP-MAS

Masimo Ear Probe  
1895-LNCS-TC1



### SP02MP-MAS

SpO2 Multisite Reusable Sensor Masimo  
2258-LNCS-Y1



### SP02IC-MASM20-04

SpO2 Interface Cable for Masimo  
4252-RED-LNC-M20

## Viamed


**SP02DC-VE6500** SpO2 Lingual Ear Probe  
for use with Viamed Pulse Oximeter.



## Fukuda-Denshi

### SP02DP-FD

SpO2 Lingual Ear Probe for use with Fukuda Denshi Monitors. DEC10.

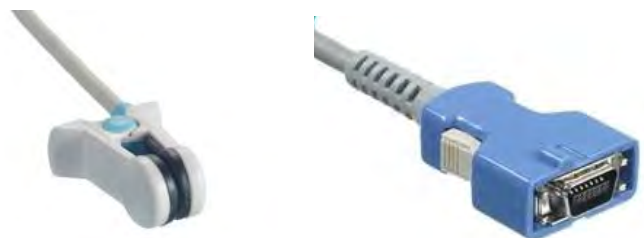

01491 527027



# Ambu

- Highly conductive wet gel - ensures reliable traces
- Superior adhesion - prevents false alarms
- Comfortable foam backing - is easy to use
- Unique offset connector - prevents artifacts from disrupting readouts
- High quality Ag/AgCl sensor - ensures excellent trace quality



**ECGP5M (N-00-S)**  
Ambu® BlueSensor N  
Pack of 25



**ECGPMED (M-00-S)**  
Ambu® BlueSensor M  
Pack of 50



**ECGPLA (L-00-S)**  
Ambu® BlueSensor L  
Pack of 25



**ECGPXLA (VL-00-S)**  
Ambu® BlueSensor VL  
Pack of 25

**ECGP/R-00-S**  
Ambu® BlueSensor R  
Pack of 25



**ECGP4MM**  
Ambu® BlueSensor M  
Pack of 50



The Ambu WhiteSensor 4570M features a wet gel with good adhesion to ensure a good signal quality during paediatric short-term monitoring ECG. Thanks to the flexible foam backing material and the ergonomic shape, the electrode ensures ease of use and comfort during application. It is radiolucent and MR Conditional.

**ECGMRI-4570M**  
Ambu® White Sensor  
Pack of 30



**ECGMRI (40713)**  
Ambu® White Sensor  
Pack of 30



**MRI SAFE**

**MRI SAFE**

01491 527027



Skintact FS-50 pressure electrodes. Designed for rest and exercise ECG.

"Aqua wet" adhesive gel for better grip.

**ECGFS50**

Skintact Electrodes 50mm dia



Defibrillator Pads



Code	Description
<b>EGCLPE</b>	Physio Control LIFEPAK EDGE System Electrodes with QUIK-COMBO connector system 1 x set of electrodes (11996-000093)
<b>EGCLPE-A</b>	Physio Control LIFEPAK EDGE System Electrodes with QUIK-COMBO connector system (adult) 1 x set of electrode (1196-000017)

Compatible with:  
LIFEPAK 12, LIFEPAK 15, LIFEPAK 20



Code	Description
<b>ECG-MR62</b>	Mindray BeneHeart MR62 Adult/Paediatric Electrode Pads
<b>EGC-MR63</b>	Mindray BeneHeart MR63 Adult/Paediatric Electrode Pads



01491 527027

## ECG MRI Non-Magnetic Quadtrode

Code	Description
<b>ECGPQ-A</b>	Adult/Paediatric Quadtrodes Pack of 25
<b>ECGPQ-N</b>	Neonate Quadtrodes Pack of 25

Four electrodes in a single, easy to use disposable patch

Designed for general adult MRI use with a large electrode contact area to provide optimal ECG performance. This unique design incorporates four electrodes into a single, easy to use disposable patch. Easy to apply, single disposable pad.

Designed Specifically for MR Use  
 Disposable Line for infection control  
 Reusable Line for reduced operating costs

Provides optimal performance and safety  
 Reduces patient preparation time



## ECG Clips

Code	Description
<b>ECGCC</b>	Comfy Clips. Pack of 4
<b>ECGCIMEC</b>	Veterinary Alligator Clip for IMEC. Pack of 5



Comfy Clips



IMEC Clips

## ECG Adapters

Code	Description
<b>ECGBC4</b>	4mm Banana to Clip, Set of 10
<b>ECGBC3</b>	3mm Banana to Clip, Set of 10
<b>ECGBSN</b>	Banana Clip to Snap. Set of 10
<b>ECGBPA</b>	Banana Plug Adapters. Pack of 4



ECGBC4 & ECGBC3



ECGBSN



ECGBPA

## ECG Gel

Code	Description
<b>ECG-EG60G</b>	Electrode Gel Spectra 360 60g
<b>ECG-EG250G</b>	Electrode Gel Signagel 250g
<b>NUPREP</b>	Nuprep Skin Prep Gel 114g. Pack of 3



01491 527027

## Datex-Ohmeda

### ECGTC-TX-2395-I0

ECG Trunk Cable for Datex/GE (545300).  
To be used with **ECGL-LX3-90P-I0** or **ECGL-LX390S-I0**



### ECGL-LX3-90P-I0 Pinch Clamps


#### ECGL-LX3-90P-I0 Snap

ECG Lead Wire for Datex/GE (545315). To be used with **ECGTC-TX-2395-I0**



### ECG-TQ-2586-I0


ECG Trunk Cable for Datex (2017003-003)  
for use with **ECGL-LQ-90P-I0** or **ECGL-LQ-90S-I0**

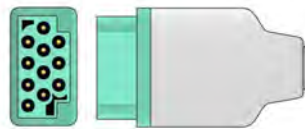

### ECGL-LQ-90P-I0 Pinch Clamps

#### ECGL-LQ-90S-I0 Snap

ECG Lead Wire for Datex (2017003-003) for  
use with **ECG-TQ-2586-I0**

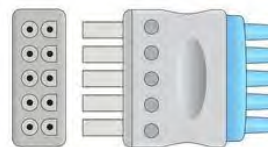

### ECGL-C2395P-I0 Pinch Clamps

#### ECGL-C2395S-I0 Snap

ECG 3 Lead Cable with Pinch Clamps for  
Datex.



### ECGL-C2386P-I0 Pinch Clamps

#### ECGL-C2386S-I0 Snap Leads

ECG 3 Lead Cable with Pinch Clamps for  
Datex.



### ECGTC-10315

ECG Trunk Cable for Datex



01491 527027

# Mindray

**ECG3L-C2312P-I0 Pinch Clamps**

**ECG3L-C2312P-I0 Snap**

ECG Single Cable for Mindray. 12 Pin.



**ECGTC-TR-2512-I0**

ECG Trunk Cable for Mindray.


To be used with ECGLW-LR3-90P0



**ECGLW-LR3-90P0 Pinch Clamps**

**ECGLW-LR3-90S0 Snap**

ECG Lead Wire for Mindray. To be used with ECGTC-TR-2512-I0

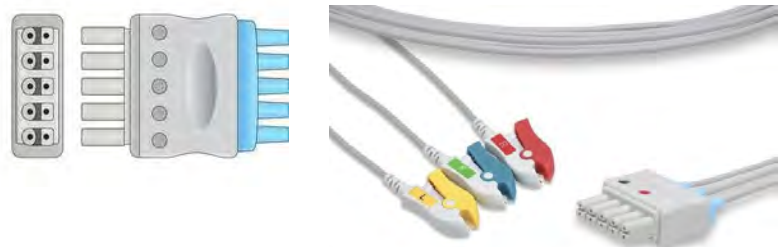

**ECGL-C2345S-I0**

ECG Single Cable for Mindray. 3 Lead. 6 Pin



01491 527027

## Fukuda Denshi

**ECGL-C2359P-I0 Pinch Clamps**

**ECGL-C2395S-I0 Snap**


ECG Direct 3 Lead Cable for Fukuda Denshi. .



**ECGL-LMB3-90P-I0 Pinch Clamps**

**ECGL-LMB3-90S-I0 Pinch Snap**

ECG Leadwire for Fukuda Denshi. 3 Leads. Use with ECG-TM-2599-I0

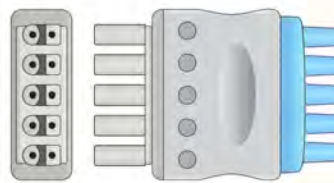

**ECG-TM2599-I0**

ECG Trunk Cable for Fukuda Denshi. Use with ECGL-LMB3-90-I0




## Nihon Kohden

**ECGL-BR903P Pinch Clamps**

**ECGL-B903S Snap**

ECG Leadwire for use with ECGTC-JC-906PA

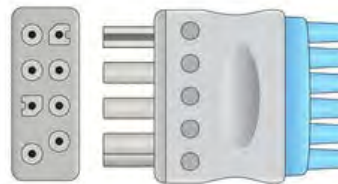

**ECG-JC-906PA**

ECG Trunk Cable for Nihon Kohden



01491 527027

## NIBP Cuffs

Available individually or in boxes of 20

Supplied with Male Luer Connectors

Larger sizes available.



Soft



Vinyl

Code	Lumen	Size	Box Qty	Box Code
<b>NIBPCDT1</b>	Double	Size 1 Soft 3-6cm	20	<b>NIBPCDT1/20</b>
<b>NIBPCDT1V</b>	Double	Size 1 Vinyl 3-6cm	20	<b>NIBPCDT1V/20</b>
<b>NIBPCST1</b>	Single	Size 2 Soft 4-8cm	20	<b>NIBPCST2/20</b>
<b>NIBPCST1V</b>	Single	Size 2 Vinyl 4-8cm	20	<b>NIBPCST2V/20</b>
<b>NIBPCDT2</b>	Double	Size 2 Soft 4-8cm	20	<b>NIBPCDT2/20</b>
<b>NIBPCDT2V</b>	Double	Size 2 Vinyl 4-8cm	20	<b>NIBPCDT2V/20</b>
<b>NIBPCST2</b>	Single	Size 2 Soft 4-8cm	20	<b>NIBPCST2/20</b>
<b>NIBPCST2V</b>	Single	Size 2 Vinyl 4-8cm	20	<b>NIBPCST2V/20</b>
<b>NIBPCDT3</b>	Double	Size 3 Soft 6-11cm	20	<b>NIBPCDT3/20</b>
<b>NIBPCDT3V</b>	Double	Size 3 Vinyl 6-11cm	20	<b>NIBPCDT3V/20</b>
<b>NIBPCST4</b>	Single	Size 4 Soft 7-13cm	20	<b>NIBPCST4/20</b>
<b>NIBPCST4V</b>	Single	Size 4 Vinyl 7-13cm	20	<b>NIBPCST4V/20</b>
<b>NIBPCDT4</b>	Double	Size 4 Soft 7-13cm	20	<b>NIBPCDT4/20</b>
<b>NIBPCDT4V</b>	Double	Size 4 Vinyl 7-13cm	20	<b>NIBPCDT4V/20</b>
<b>NIBPCST5</b>	Single	Size 5 Soft 8-13cm	20	<b>NIBPCST5/20</b>
<b>NIBPCST5V</b>	Single	Size 5 Vinyl 8-13cm	20	<b>NIBPCST5V/20</b>
<b>NIBPCDT5</b>	Double	Size 5 Soft 8-13cm	20	<b>NIBPCDT5/20</b>
<b>NIBPCDT5V</b>	Double	Size 5 Vinyl 8-13cm	20	<b>NIBPCDT5V/20</b>
<b>NIBPCST6</b>	Single	Size 6 Soft 12-19cm	20	<b>NIBPCDT6/20</b>
<b>NIBPCST6V</b>	Single	Size 6 Vinyl 12-19cm	20	<b>NIBPCDT6V/20</b>
<b>NIBPCDT6</b>	Double	Size 6 Soft 12-19cm	20	<b>NIBPCDT6/20</b>
<b>NIBPCDT6V</b>	Double	Size 6 Vinyl 12-19cm	20	<b>NIBPCDT6V/20</b>



Adapts to 1-tube or 2-tube BP monitors

Code	Size
<b>NIBPCDR1</b>	Size 1 3-6cm
<b>NIBPCDR2</b>	Size 2 4-8cm
<b>NIBPCDR3</b>	Size 3 6-11cm
<b>NIBPCDR4</b>	Size 4 7-13cm
<b>NIBPCDR5</b>	Size 5 8-13cm



01491 527027

## Philips Compatible



Code	Size	Box Qty	Box Code
NIBPCPH1	Size 1 3-6cm	24	NIBPCPH1/24
NIBPCPH2	Size 2 4-8cm	24	NIBPCPH2/24
NIBPCPH3	Size 3 6-11cm	24	NIBPCPH3/24
NIBPCPH4	Size 4 7-13cm	24	NIBPCPH4/24
NIBPCPH5	Size 5 8-13cm	24	NIBPCPH5/24



Fits onto these connectors

## NIBP Accessories

Code	Description
NIBPACTT	Blood Pressure Airhose Tubing 3.2mmID/6.4mm OD. 1 metre
NIBPLSC	NIBP Male Slip Connectors
NIBPLSCFM	NIBP Luer Slip Connector Female, Metal
NIBPHP-MR	NIBP Hose Pin for Mindray

We are able to source many different NIBP Cuff Connectors, so if you need one that is not listed here, please contact us on 01491 527027 or sales@acevetssupplies.co.uk.



Male Slip Luer



Female Slip Luer



Hose Pin for Mindray

## Aneroid Sphygmomanometer Palm Model

- Bright, easy to read dial surface housed in tough ABS plastic
- One handed control for easy air release with either hand
- Precision mechanism ensures smooth performance & immediate response
- Compact enough to carry in zip-up case supplied
- Measurement range: 0-300mmHg
- Accurate to +3mmHg
- 1 year warranty



Code	Description
AS-BK2012	Aneroid Sphygmomanometer

01491 527027



## GE/Datex-Ohmeda

### **NIBPH2017009-001**

1 Piece NIBP Hose for GE Marquette.  
Neonatal 2.4M



### **NIBPH2017009-003**

1 Piece NIBP Hose for GE Marquette.  
Neonatal 3.6M



### **NIBPH877514**

1 Piece NIBP Hose for GE Marquette.  
Neonatal 2.4M



## Mindray

### **NIBPH6200-30-09688**

1 Piece NIBP Hose Mindray/  
Datascop



### **NIBPH-ASN-15-200**

1 Piece NIBP Hose Mindray/  
Datascop



## Fukuda-Denshi

### **NIBPH-AS-04-200**

1 Piece NIBP Hose Fukuda Denshi and  
Mindray/Datascop. 2.5M



Please ask us for any other NIBP Cables you may need.

Also for extra long MRI cables.

**01491 527027**

## IBP Transducer Cables

We can source many manufacturer's IBP cables some of which are detailed below.

Please advise which Transducers you use and the make of your Multi Parameter Monitor.

Contact our sales office on 01491 527027 or email [sales@acevetssupplies.co.uk](mailto:sales@acevetssupplies.co.uk)



ARGON



B. BRAUN



BD



EDWARDS




MEDEX ABBOTT



MEDEX LOGICAL



UTAH

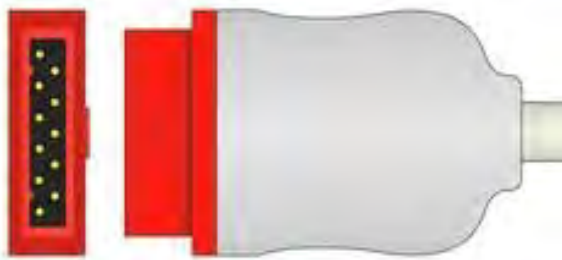

### GE—MARQUETTE

**IC-MQ-ED0** Edwards

**IC-MQ-BD0** BD

**IC-MQ-MX0** Medex Abbott

**IC-MQ-MX10** Medex Logical



### IBPDAC-GE Dual Adapter



### DATEX OHMEDA

**IC-DX-ED0** Edwards

**IC-DX-BD0** BD

**IC-DX-MX0** Medex Abbott



### FUKUDA-DENSHI



**IC-FD-ED0** Edwards

**IC-FD-BD0** BD

**IC-FD-MX0** Medex Abbott

01491 527027

# IBP Transducer Cables

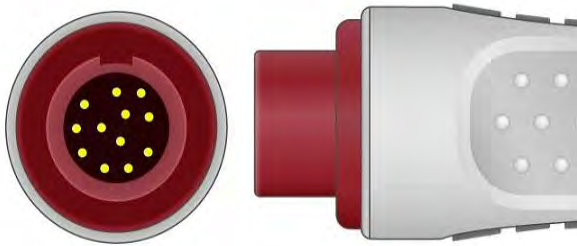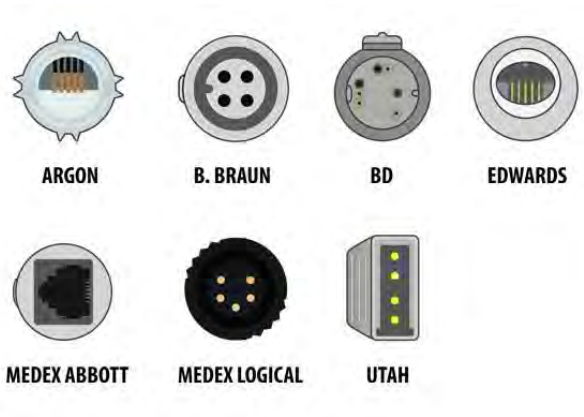

**MINDRAY-DATASCOPE (12 PIN)**

- IC-MR-ED0**      Edwards
- IC-MR-BD0**      BD
- IC-MR-MX0**      Medex Abbott 001C-30-70759



**IBPDAC-MR12P** Mindray 12 Pin Dual Adapter



**MINDRAY-DATASCOPE (6 PIN)**

- KA-6P-ED**      Edwards

01491 527027


## IBP Transducer Kits



**EDWARDS**

**IP-ED-300/20 Box of 20**


**IP-ED-300 each**

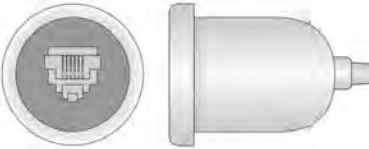

**MEDEX ABBOTT**

**IP-MX-300/20 Box of 20**


**IP-MX-300 each**

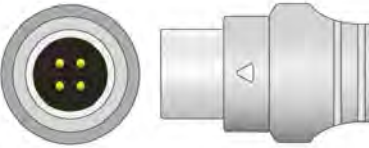

**B. BRAUN**

**IP-BB-300/20 Box of 20**

**IP-BB-300 each**

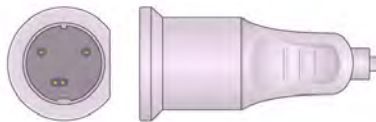

**BD**

**IP-BD-300/20 Box of 20**

**IP-BD-300 each**

**SMITHS MEDICAL LOGICAL SINGLE PRESSURE MONITORING KIT, CLEAR, 210CM, 180CM PLUS 30CM EXTENSION.**

PMK-MX9605/20 BOX OF 20

PMK-MX905 each

**PRESSURE TRANSDUCER PLATE**

PT-MX9601P1CZ each



01491 527027

## Vetronic MiniTemp

Code	Description
<b>MINITEMP</b>	Vetronic Minitemp Unit, includes choice of probe
<b>MT400RP</b>	Rectal Probe for 400 Series
<b>MT400SP</b>	Surface Probe for 400 Series
<b>MT415RP</b>	Rigid Probe for 415 Series
<b>MT415FP</b>	Flexible Probe for 415 Series

The MiniTemp is a small, hand-held battery operated temperature monitor designed specifically for use with animals. The monitor supports two common ranges of temperature sensor – the industry-standard YSI 400-compatible range and the 8415 miniature range. The miniature sensors measure just 1.2mm in diameter and are suitable for use with small mammals and exotic species. The standard temperature sensors measure 4mm in diameter and are suitable for all larger species. The introduction of the concept of a comfort zone for temperature monitoring is designed to enable the user to see at a glance if the animal is within the comfortable range of body temperatures. The comfort zone is designed for use with mammals and has a lower limit of 36.6°C (97.9°F) and an upper limit of 39.6°C (103.2°F). This relates to one degree above and two degrees below normal body temperature. When the temperature being monitored is within the comfort zone the LED at the top of the unit flashes green. For temperatures outside the comfort zone, the LED flashes red. The unit is housed in a robust rubber boot resistant to drops and knocks. The boot also has a built-in tilt stand for improved viewing angles.

### Mechanical

- Weight 125g (225g with -fitted rubber boot)
- Dimensions 86mm x 145mm x 34mm
- Construction Rugged ABS plastic case housed in a rubber boot

### Electrical

- Screen size 25mm x 18mm
- Screen type 128 x 64 pixels, transfective LCD
- Electrical supply Alkaline or rechargeable 9v PP3 battery only.
- Battery life 200 hours run time, six months standby

### Temperature Sensors\*

Temperature range 5°C to 45°C

Accuracy +/- 0.1°C within specified range. Accuracy outside of this range is unspecified

Standard 400 Series:

Rectal type 2m x 4mm, terminated in 6.35mm standard jack plug

Surface type 2m x 8mm disc, terminated in 6.35mm standard jack plug

Miniature Series (8514):

Rigid probe 1.2mm tip housed in a 2mm x 80mm rigid metal probe

Flexible probe 1.2mm tip housed in a 2mm x 270mm flexible sheath, terminated with a 3.5mm jack plug.

\*MiniTemp can be fitted with a module to support either the standard sensors, or the miniature sensors, but not both. Requirements must be specified when ordering.



01491 527027

## Veterinary Thermometer

- 10 second measurement time
- Auto shut off after 10 minutes
- 32.0°C to 42.9°C measurement range
- Clear digital changeable display ( °C / °F )
- Waterproof



Easy to use battery operated electronic thermometer with chart giving Hobday's specification (the normal body temperature for a variety of animals).

The temperature range of this thermometer is 32 to +42°C and it displays the temperature in 0.1°C increments. It has an accuracy of +/-0.1°C between 35 and 39°C, otherwise +/-0.2°C.

It uses a single LR41 or equivalent battery.

Code	Description
<b>THERMVET</b>	Veterinary Thermometer

## Flexible Tip Thermometer

- Flexible tip
- Auto shut off
- 1.5v button battery
- Last memory reading
- Bleeper function
- Size: 134 x 22 x 14mm



Code	Description
<b>THERMFLEX</b>	Flexible Tip Thermometer
<b>THERMPC100</b>	Thermometer Probe Cover. Pack of 100



## Infrared Tympanic Ear Thermometer

- Quick and Accurate measuring time <1 second.
- Long life: over 4000 readings with 2x AAA batteries.
- Automatic Power off.
- Display range: 34°C – 44°C.
- Stores last 10 readings.
- Low battery warning.



Code	Description
<b>THERMEAR</b>	Ear Thermometer with Case
<b>THERMEPC20</b>	Ear Thermometer Probe Covers. Pack/ 20.


01491 527027

# Temperature Probes


Code	Description
<b>TPMR401</b>	Temperature Probe with Jack Plug Large
<b>TPMR402</b>	Temperature Probe with Jack Plug Small



Code	Description
<b>TPMR401B</b>	Temperature Probe Large for Mindray Beneview
<b>TPMR402B</b>	Temperature Probe Small for Mindray Beneview

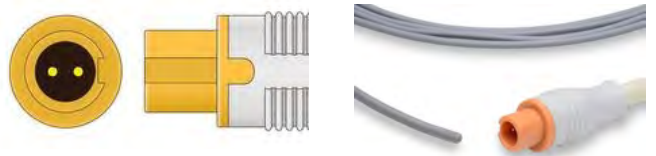


Code	Description
<b>TP-DMQ-AD0</b>	Temperature Probe Large for Datex/GE 11 Pin Flat

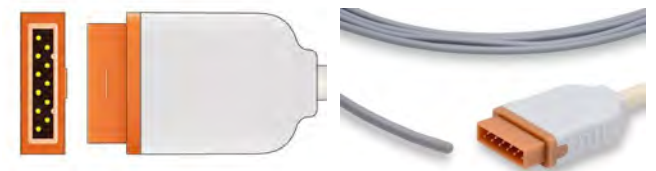

Code	Description
<b>TS-D2252-AS0</b>	YSI Compatible Reusable Temperature Probe
<b>TS-T004-160</b>	Compatible Disposable, Adhesive Temperature Probe Covers. Pack of 16.



Code	Description
<b>TPA-DMQ-AG0</b>	Temperature Probe Adapter Cable for Datex/GE 11 Pin Flat. To allow use Jack Plug Probe.

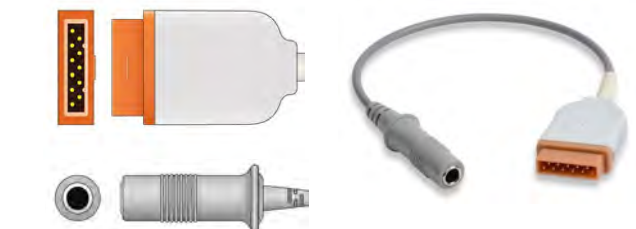

01491 527027

## Gas Sample Lines

Available individually or in box quantities

Code	Ends	Length	Box Code	Box Qty
<b>GSLMF3</b>	Male/Female	3 metres	<b>GDLMF3/50</b>	50
<b>GSLMM2GE</b>	Male/Male. GE Healthcare	2 metres	<b>GSLMM2GE/10</b>	10
<b>GSLMM3</b>	Male/Male	3 metres	<b>GSLMM3/50</b>	50
<b>GSLMM3GE</b>	Male/Male. GE Healthcare	3 metres	<b>GSLMM3GE/10</b>	10
<b>GSLMM6GE</b>	Male/Male. GE Healthcare.	6 metres	<b>GSLMM6GE/5</b>	5
<b>GSLMM9</b>	Male/Male.	9 metres		
<b>LUERFF</b>	Female to Female Luer Connector			
<b>SSA2215</b>	Side Stream Adapter 22/15mm Female with Sample Port. Dead space 6.6ml.			



LUERFF Female to Female Luer Connector (use on GSLMF3 to achieve 2 female ends)



SSA2215

Single Code	Description	Length	Box Code	Box Qty
<b>GSLMS-Y</b>	Gas Sampling Line with T Connector, Microstream, Yellow. Long Term—up to 72 hours.	2.1 metres	<b>GSLMS-Y/25</b>	25
<b>GSLMS-O</b>	Gas Sampling Line with T Connector, Microstream, Orange. Short Term—up to 8 hours	2.1 metres	<b>GSLMS-O/25</b>	25

Designed to connect to low flow (50ml/min) Microstream™ enabled monitors and provide a crisp, accurate waveform

- Full range of sampling options including single nostril, nasal cannula and intubated airway adapter to cover a wide range of patients and applications
- Colour coded packaging to make selecting the right item easy
- 0.2-micron filter included on all lines to protect the monitor from contamination with streamlined design to avoid patient discomfort



01491 527027



# Water Traps and Filters



**MAWT-8208**  
Surgivet Water Trap 8208.



**MASRP15**  
Monitoring Line Hydrophobic Filter with Male and Female Luer.



**MAWTDOBK**  
D'Fend Water Trap, Datex-Ohmeda. Black 876446-HEL



**MAWTDOMI**  
D'Fend Water Trap, Datex-Ohmeda. Mini. 8002174



**MAWT - M1182629**  
D-Fend Pro Water Trap. Datex-Ohmeda Grey.



**MAWT-M1200227**  
D-Fend Pro Water Trap. Datex-Ohmeda Green.



**MAWTMR**  
Mindray Adult/Paediatric Water Trap. 60-13100-00



**MAWTM02C**  
Mindray Adult/Paediatric Water Trap.



**MAWTMRN**  
Mindray Neonatal Water Trap. 100-000080-00



**MAWT-8075**  
Surgivet Water Trap 8075.



**MAWT-3010487**  
MEDRAD Water Trap. 3010487



**CR664523**  
Water Trap Adapter Microflow Sample Line

01491 527027

## D-Lite Connectors & Tubing, GE

Code	Size
<b>MA DLITE</b>	D-Lite Sensor
<b>MA PLITE</b>	Pedi-Lite Sensor



Paediatric Sensor

Code	Size
<b>MADLST2M</b>	D-Lite Sensor Tubing 2 Metres
<b>MADLST3M</b>	D-Lite Sensor Tubing 3 Metres



## Intersurgical Spirometry Sets

Code	Size
<b>INT2740003</b>	Adult Spirometry Set 3 metres
<b>INT2740006</b>	Adult Spirometry Set 6 metres
<b>INT2741003</b>	Paediatric Spirometry Set 3 metres
<b>INT2741006</b>	Paediatric Spirometry Set 6 metres



Adult Spirometry Set

Code	Size
<b>INT2742000</b>	Adult Spirometry Sensor
<b>INT2743000</b>	Paediatric Spirometry Sensor



Paediatric Spirometry Set



Adult Sensor



Paediatric Sensor

## Manometer Lines

Code	Size
<b>ML210MFL</b>	Manometer Line 210mm Male/Female Luer Connectors



01491 527027

## Vetronic ICON Oxygen Analyser

Icon is a new in-circuit oxygen monitor designed for use in veterinary operating theatres and intensive care wards.

The Icon unit can be placed on the Fresh Gas Outlet of an anaesthetic machine or within the rebreathing limbs of a circle system, wherever the oxygen concentration needs to be known.

The Icon unit uses a single AA battery that has a lifetime of approximately 1 year when used for 8 hours a day every day. The sensing element is a longlife galvanic sensor that has a working life of approximately 2 years.



Code	Description
ICON	ICON Oxygen Analyser

## Datex-Ohmeda Oxygen Cell

Code	Description
M-10	M-10 Oxygen Cell 6050-0004-110



Other types available on request.  
Email  
sales@acevetsupplies.co.uk or  
call 01491 527027.

## Oesophageal Stethoscopes

Code	Size	Code	Size
OS12SH	12fr x 500mm	OS12LO	12fr x 1100mm
OS18SH	18fr x 500mm	OS18LO	18fr x 1100mm
OS24SH	24fr x 500mm	OS24LO	24fr x 1100mm



### Ear Piece for Oesophageal Stethoscope

Code	Size
OSSTETH	Ear Piece only
SET2	Set of 2 Soft Littman Eartips



## FOB Watches

- New improved design
- Japanese Quartz movement
- Pin clasp closure
- Made from silicone based material
- Ideal for infection control
- Contains 1 x watch battery (SR626SW)
- 2 year battery guarantee



Code	Size
FOBWAT	Please specify colour
FOBWAT14	Pack of 14 (2 of each colour)

01491 527027

## Suction Receiver Jar, Autoclavable

Autoclavable polysulphone collection jar, complete with cap assembly and overflow safety valve.

Code	Description
<b>SRJ15</b>	1500ml
<b>SRJ20</b>	2000ml
<b>SRJRC15</b>	1500ml with Rail Clamp
<b>SRJRC20</b>	2000ml with Rail Clamp
<b>SRJVW15</b>	1500ml with 'V' Wedge
<b>SRJVW20</b>	2000ml with 'V' Wedge



'V' Wedge

Rail Clamp

## Suction Catheter, PVC. Pack of 10

Code	Size	Code	Size
<b>SC0610</b>	6fr	<b>SC1410</b>	14fr
<b>SC0810</b>	8fr	<b>SC1610</b>	16fr
<b>SC1010</b>	10fr	<b>SC1810</b>	18fr
<b>SC1210</b>	12fr		



- Manufactured from transparent medical grade PVC allowing clinicians to visualise any blockages quickly.
- Atraumatic tip with double offset eye to provide relief if the distal hole becomes occluded.
- Smooth surface for easy insertion.
- Different connector ends available for compatibility with a range of suction units and aspirators.
- Colour coded according to size to make visual identification easy.
- Supplied sterile in individual easy access sachets.

## Suction Tubing

Code	Size
<b>ST6518</b>	6.5mm ID x 1.8M Length
<b>ST6530</b>	6.5mm ID x 3.0M Length

- Transparent tubing allows blockages to be visualised quickly.
- Thick walls prevent the tube from collapsing under high negative pressure.
- Flexible universal connectors provide a secure connection to the handle and source.
- A range of lengths available.
- Supplied sterile in individual easy access sachets.



01491 527027

## Suction Yankauer Handle

Available individually or in boxes of 20

Code	Description	Box Qty	Box Code
<b>SYVC431001</b>	Suction Yankauer Handle Vented Crown Tip	20	<b>SYVC431001/20</b>
<b>SYVP431004</b>	Suction Yankauer Handle Vented Plain Tip	20	<b>SYVP431004/20</b>



- Transparent suction handles for clear visualisation during suction. Similar in shape to traditional metal Yankauer suckers.
- Available with plain or crown tip, with or without vents to control suction.
- Also available packed with suction tubing and handles for maximum convenience.
- Supplied sterile in individual easy access sachets.

## Manual Suction Pump

Code	Size
<b>SPMCOM</b>	Suction Pump Complete
<b>SPMRDJ</b>	Replacement Jar with Adult & Infant Catheters

- Hand powered, can be used anywhere
- Portable, compact and lightweight
- Easy to use with one hand
- Auto stop when full



## Suction Controller

The suction controller enables the vacuum to be set at any level from 0 to full pipeline vacuum. The vacuum gauge is marked from 0 to -100 kPa and 0 to -750 mmHg.

Code	Size
<b>SC0100</b>	Suction Controller with Direct BS Probe



## Eschmann Suction Filter

Code	Size
<b>SF-VP25</b>	Suction Filter for Eschmann VP25, VP35, VP45. Box of 10.



01491 527027

## Waste Gas Canisters

Code	Description
<b>WAG001</b>	Waste Gas Canister Single
<b>WAG006</b>	Waste Gas Canister Box of 6

- Contains activated carbon
- Absorbs 200g of organic waste anaesthetic agent
- Not suitable for use with N<sub>2</sub>O
- Approximate maximum constant use—20 hours
- Discard when weight reaches 1400g
- Dispose by controlled incineration



## Connectors & Tubing



INT1971



INT1970



INT1971001



INT1961



GE2105488-016



Code	Size
<b>INT1961</b>	22F-22M (Connects 22mm tubing to the top of the WAG Canister)
<b>INT1970</b>	22M-30M
<b>INT1971</b>	22M-30F
<b>INT1971001</b>	30M-30F
<b>GE2105488-016</b>	Gas Return Adapter 22mm M—22mm F



Code	Size
<b>GEL-M1031274</b>	Exhaust Line for gas return and scavenging

Code	Size
<b>INT1574</b>	22mm x 400mm Cuffed Length

## APL Valve

Code	Description
<b>INT1972000</b>	Anaesthesia gas scavenging pressure relief valve 30F/30F with 4.2M 30mm Flextube
<b>INT1972001</b>	Anaesthesia gas scavenging pressure relief valve 30F/22M with 22mm Flex-tube, ≥ 4.8m

This in-line pressure relief valve is connected to the outlet port of the breathing system APL valve. The purpose of the valve is to allow a pressure release in the case of the AGSS transfer tube becoming occluded



01491 527027

## MEC Receiving System

Conforms to ISO 8835-3

- MRI compatible
- Wall, rail or machine mounted
- Weighs 1 Kg, size 80mm X 390mm
- Supplied with Disposal and Transfer hoses
- Both Disposal Hose and Patient Inlet port swivel about 360° to ensure hose positioning and prevent occlusion
- 3m, 4m or 5m Disposal Hose
- 2.4m Transfer Hose
- Complete with Pressure Relief Valve.



Code	Size
12800000	Receiving System inc. 5m Disposable Hose and 2.4m Transfer c/w Pressure Relief Valve
12800013	Replacement Hose Assembly

## Purair 80 AGSS

Code	Size
PURAIR80	Purair 80 Gas Receiving System



Designed as a stand alone scavenging unit for use with anaesthetic work stations that have an integrated patient air break

- Cost effective method of anaesthetic gas scavenging without the need for remote pipework
- Simple to install and factory set flow rates, so no commissioning is required
- Extraction capacity is 80 litres per minute to current European flow rates
- Conform to ISO 8835-3 and is used for the extraction of waste and excess gases
- Occlusion and fault alarm.

01491 527027

## Key Fill Bottle Adapters

Code	Description
<b>KFBA-ISO</b>	Key Filler for Isoflurane
<b>KFBA-SEV</b>	Key Filler for Sevoflurane



## Pour Fill Bottle Adapters

Code	Description
<b>PFBA-ISO</b>	Pour Filler for Isoflurane
<b>PFBA-SEV</b>	Pour Filler for Sevoflurane



## Easy-Fil Bottle Adapters

Code	Description
<b>EFBA-ISO</b>	Key Filler for TEC 7, TEC 8 and Aladdin 2 cassette -Isoflurane
<b>EFBA-SEV</b>	Key Filler for TEC 7, TEC 8 and Aladdin 2 cassette - Sevoflurane



## 'O' Rings

Code	Description
<b>ORINGVAP</b>	Vaporiser O Ring



Please call or email our Sales Department on sales@acevetssupplies.co.uk or 01491 527027 for vaporiser prices.

01491 527027






Code	Description
<b>INT21223M</b>	3M Long Ayres T Piece with APL Valve & 0.5 litre bag
<b>INT21226M</b>	6M Long Ayres T Piece with APL Valve & 0.5 litre bag




Code	Description
<b>CIRCSBC1530</b>	Smooth Bore Y Circuit 15mm x 3m with 22F end connections
<b>CIRCSBC1560</b>	Smooth Bore Y Circuit 15mm x 6m with 22F end connections (Connections at 3m)

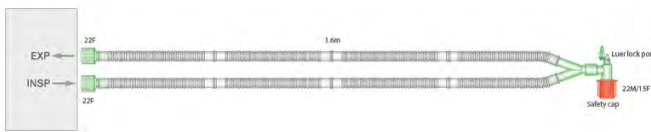

Code	Description
<b>CIR2232</b>	Circle Circuit 22mm x 3.2m
<b>CIR1532</b>	Circle Circuit 15mm x 3.2m
<b>CIR2260</b>	Circle Circuit 22mm x 6m
<b>CIR1560</b>	Circle Circuit 15mm x 6m

These circuits can be made to any length required. Call 01491 527027 or Email sales@acevetssupplies.co.uk for details

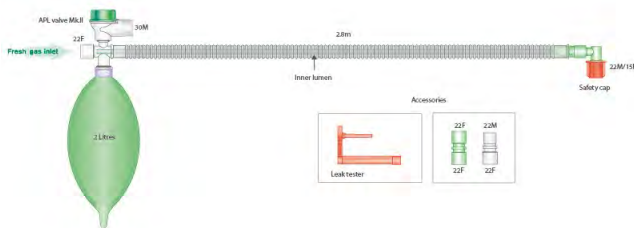

Code	Description
<b>BAIN54</b>	Bain 22mm x 5.4m (Intersurgical)



Code	Size
<b>GSLMM6GE</b>	Male/Male Luer Connectors 6M
<b>GSLMM9</b>	Male/Male Luer Connectors 9M



Code	Size
<b>INT2740003</b>	Adult Spirometry Set 3 metres
<b>INT2740006</b>	Adult Spirometry Set 6 metres
<b>INT2741003</b>	Paediatric Spirometry Set 3 metres
<b>INT2741006</b>	Paediatric Spirometry Set 6 metres

01491 527027



Code	Description
<b>ECGPQ-A</b>	Adult/Paediatric Quadrodes Pack of 25
<b>ECGPQ-N</b>	Neonate Quadrodes Pack of 25

Full details on Page 73



**ECGMRI-4570M**  
Ambu® White Sensor  
Pack of 30

**MRI SAFE**



**ECGMRI**  
Ambu® White Sensor  
Pack of 30

**MRI SAFE**

## MRI Safe Thermal Mattress

Designed exclusively for use in MRI  
 No MRI Signal Loss  
 Clinically Tested to Maintain Body Temperature  
 Available in 3 Sizes

- Small 510mm x 710mm
- Medium 610mm x 860mm
- Large 710mm x 1015mm

Two Layer Internal Design  
 Insulative conductive thermal barrier  
 Proprietary MRI Safe heat loss barrier



Full details on Page 49

## MRI Compatible Disposable Laryngoscopes




Full details on Page 35

01491 527027

## Syringe Labels

Code	Description
SLWMT	Wall Mounted Holder for 10 boxes of Syringe Labels
SLTT	Tidy Tray for 10 boxes of Syringe Labels


- Any colour and any wording available.
- Over 200 types held in stock
- No minimum order
- Non standard labels subject to a £25.00 plate charge

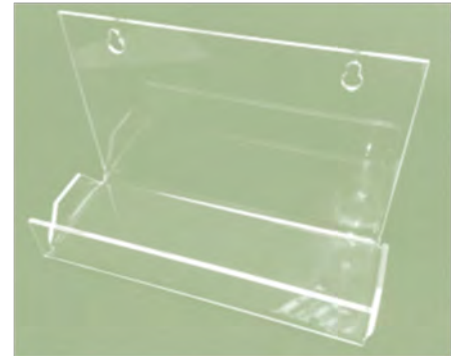

Please see back pages for a full list of our syringe labels.

## Dose Labels

Code	Description
LAB-DOSE	Clear Vinyl Dose Label. Roll of 500.

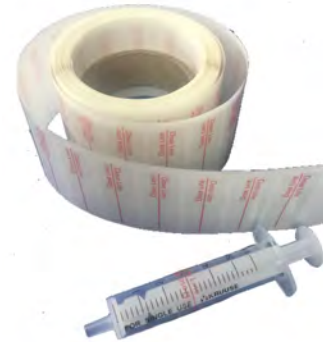

## Monitoring Labels

Code	Description
LAB-3029	Do Not Pull Me. Strip of 20
LAB-3022	Do Not Throw Away. Strip of 20



## Information Labels



Code	Description
LAB-INFE	Infection Risk, Green Label. Roll of 1,000
LAB-DNR	D.N.R. Red Label. Roll of 1,000



01491 527027

## Other Labels

Code	Description
LAB-FEED	Feeding Tube Label. Roll of 1,000. 12cm x 2.5cm
LAB-DRA100	Drugs Added Label. Pack of 100. 8cm x 4cm
LAB-G-FLUID	Fluid Type Labels, Green. Roll of 1,000. 6.2cm x 3.8cm
LAB-G-SYRI	Drug Detail Labels, Green . Roll of 1,000. 6.2cm x 3.8cm
LAB-FLUIDS	Fluids Signed Label, White. Roll of 1,000. 5cm x 3.8cm

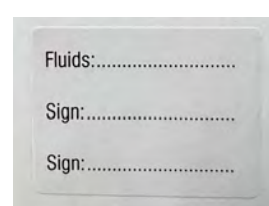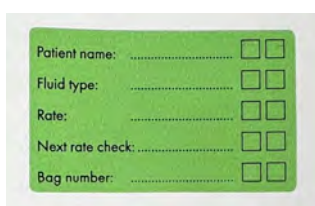

## Syringe Labels

SL-O-ACEP	ACEPROMAZINE, Orange	SL-W-ASPA	ASPARAGINASE, White
SL-S-ACEP	ACEPROMAZINE, Salmon	SL-OS-ATIP	ATIPAMEZOLE, Orange Stripe
SL-O-ACP	ACP, Orange	SL-S-ATIP	ATIPAMEZOLE, Salmon
SL-Y-ACP	ACP, Yellow	SL-W-ATIP	ATIPAMEZOLE, White
SL-G-ACPB	ACP/BUPRENORPHINE, Green	SL-WR-ATIP	ATIPAMEZOLE, White/Red Writing
SL-G-ACBU	ACP/BUTORPHANOL, Green	SL-P-ATRA	ATRACURIUM (KEEP IN FRIDGE) Pink
SL-W-ACMM	ACP/MEDETOMIDINE/METHADONE, White	SL-BP-ATRA	ATRACURIUM (KEEP IN FRIDGE) B.Pink
SL-G-ACPM	ACP/METHADONE, Green	SL-R-ATRA	ATRACURIUM, Red
SL-S-ACPM	ACP/METHADONE, Salmon	SL-G-ATRO	ATROPINE, Green
SL-V-ACPM	ACP/METHADONE, Violet	SL-W-AUGM	AUGMENTIN, White
SL-Y-ACPT	ACP/TORBUGESIC, Yellow	SL-Y-BEAD	BEADS, Yellow
SL-Y-ACPV	ACP/VETERGESIC, Yellow	SL-B-BLAN	BLANK, Blue
SL-R-ADRE	ADRENALINE, Red	SL-G-BLAN	BLANK, Green
SL-V-ADRE	ADRENALINE, Violet	SL-O-BLAN	BLANK, Orange
SL-R-ALFA	ALFAXAN, Red	SL-R-BLAN	BLANK, Red
SL-Y-ALFA	ALFAXAN, Yellow	SL-V-BLAN	BLANK, Violet
SL-B-ALFE	ALFENTANIL, Blue	SL-W-BLAN	BLANK, White
SL-Y-ALON	ALPHAXALONE, Yellow	SL-Y-BLAN	BLANK, Yellow
SL-W-AMOX	AMOXICILLIN, White	SL-W-BLEO	BLEOYCIN, White
SL-WR-ANTIB	ANTIBIOTIC, White (Red Writing)	SL-G-BU25	BUPIVACAINE 0.25%, Grey
SL-W-ANTI	ANTIBIOTIC, White	SL-G-BU50	BUPIVACAINE 0.5%, Grey
SL-O-ADAN	ANTISEDAN, Orange	SL-G-BUPI	BUPIVACAINE, Grey
SL-WR-ADAN	ANTISEDAN, White (Red Writing)	SL-B-BUPR	BUPRENORPHINE, Blue
SL-WR-ARTE	ARTERIAL, White (Red Writing)	SL-W-BUSC	BUSCOPAN, White

01491 527027

SL-S-BUDE	BUTORPH/DEXMED, Salmon	SL-G-DUPH	DUPH LA, Green
SL-B-BUTO	BUTORPHANOL, Blue	SL-RS-EDRO	EDROPHONIUM, Red Stripe
SL-G-BUTO	BUTORPHANOL, Green	SL-V-EPHE	EPHEDRINE, Violet
SL-V-CALC	CALCIUM GLUCONATE, Violet	SL-O-EPID	EPIDURAL, Orange
SL-R-CALC	CALCIUM, Red	SL-W-EPIR	EPIRUBICIN, White
SL-O-CAPR	CAPRIEVE, Orange	SL-VS-ESMO	ESMOLOL, Violet Stripe
SL-W-CARB	CARBOPLATIN, White	SL-Y-ETOM	ETOMIDATE, Yellow
SL-O-CARP	CARPROFEN, Orange	SL-RS-EUTH	EUTHANASIA, Red Stripe
SL-W-CARP	CARPROFEN, White	SL-R-EXPI	EXPIRES, Red
SL-W-CEFU	CEFUROXIME, White	SL-B-FENT	FENTANYL, Blue
SL-WR-CEFU	CEFUROXIME. White (Red Writing)	SL-OS-FLUM	FLUMAZENIL, Orange Stripe
SL-G-CEPH	CEPHALEXIN, Green	SL-W-FLUO	FLUOROURACIL, White
SL-S-CERE	CERENIA, Salmon	SL-O-FORM	FORMAL SALINE, Orange
SL-W-CHLO	CHLORPHENAMINE, White	SL-W-FRUS	FRUSEMIDE, White
SL-W-CISP	CISPLATIN, White	SL-W-FURO	FUROSEMIDE, White
SL-W-CLIN	CLINDAMYCINE, White	SL-G-GADO	GADOTERIC, Green
SL-W-COAM	CO-AMOXICLAV, White	SL-O-GAND	GANDOLINIUM, Orange
SL-W-COTR	CO-TRIMOXAZOLE, White	SL-W-GLUC	GLUCOSE, White
SL-W-COLL	COLLOID, White	SL-Y-GLUE	GLUE, Yellow
SL-G-CONT	CONTRAST, Grey	SL-G-GLYC	GLYCOPYRROLATE, Green
SL-R-CONT	CONTRAST, Red	SL-RS-GLYC	GLYCOPYRROLATE/NEOSTIGMINE, Red Stripe
SL-W-CONT	CONTRAST, White	SL-W-HALF	HALF STRENGTH, White
SL-W-CRYS	CRYSALLOID, White	SL-W-HART	HARTMANN'S, White
SL-W-CYCL	CYCLOPHOSPHAMIDE, White	SL-P-HEPF	HEP FLUSH, Pink
SL-R-CYTA	CYTARABINE, Red	SL-W-HEPF	HEP FLUSH, White
SL-W-CYTA	CYTARABINE, White	SL-B-HEPA	HEPARI/ SALINE, Blue
SL-R-CYTO	CYTOTOXIC, Red	SL-W-HEPA	HEPARI/ SALINE, White
SL-YR-CYTO	CYTOTOXIC, Yellow (Red Writing)	SL-W-HEPN	HEPARIN, White
SL-W-DACT	DACTINOMYCIN, White	SL-W-HYAL	HYALURONIDASE, White
SL-S-DETO	DETOMIDINE, Salmon	SL-Y-HYPE	HYPERTONIC, Yellow
SL-W-DEXA	DEXAMETHASONE, White	SL-WR-JUGU	JUGULAR LINE. White (Red Writing)
SL-Y-DEXD	DEXDOMITOR, Yellow	SL-Y-KETA	KETAMINE, Yellow
SL-O-DEXM	DEXMEDETOMIDINE, Orange	SL-R-LASP	L-ASPARAGINASE, Red
SL-R-DEXM	DEXMEDETOMIDINE, Red	SL-V-LEVE	LEVETIRACETAM, Violet
SL-S-DEXM	DEXMEDETOMIDINE, Salmon	SL-G-LEVO	LEVOBUPIVACAINE, Grey
SL-OO-DEXM	DEXMEDETOMIDINE, White with Orange Type	SL-G-LIDO	LIDOCAINE, Grey
SL-W-DEX5	DEXTROSE 5%, White	SL-OS-LIGH	LIGHT SENSITIVE, Orange Stripe
SL-O-DIAM	DIAZEMULS, Orange	SL-G-LIGN	LIGNOCAINE, Grey
SL-O-DIAZ	DIAZEPAM, Orange	SL-G-LOCA	LOCAL, Grey
SL-G-DISC	DISCARD DATE, Green	SL-W-MAGN	MAGNESIUM SULPHATE, White
SL-R-DISC	DISCARD DATE, Red	SL-W-MANN	MANNITOL, White
SL-Y-DISC	DISCARD, Yellow	SL-G-MARC	MARCAIN, Grey
SL-P-DONO	DO NOT FLUSH, Pink	SL-S-MARO	MAROPITANT, Salmon
SL-V-DOBU	DOBUTAMINE, Violet	SL-W-MARO	MAROPITANT, White
SL-WR-DOMK	DOM/KET/TORB, White (Red Writing)	SL-Y-MARO	MAROPITANT, Yellow
SL-WB-DOM	DOM/TORBUGESIC, White with Blue	SL-W-MBKE	MED/BUT/KET, White
SL-W-DOM	DOM/TORBUGESIC, White	SL-R-MEDE	MEDETOMIDINE, Red
SL-V-DOPA	DOPAMINE, Violet	SL-S-MEDE	MEDETOMIDINE, Salmon
SL-W-DOXA	DOXAPRAM, White	SL-OO-MEDE	MEDETOMIDINE, White with Orange Type
SL-W-DOXO	DOXORUBICIN, White	SL-W-MEDE	MEDETOMIDINE, White

SL-W-MBUP	MEDETOMIDINE/BUPRENORPHINE, White	SL-O-PREM	PREMED, Orange
SL-W-MBUT	MEDETOMIDINE/BUTORPHANOL, White	SL-G-PROP	PROPOFOL PLUS, Green
SL-Y-MEDK	MEDETOMIDINE/KETAMINE, Yellow	SL-Y-PROP	PROPOFOL, Yellow
SL-W-MEME	MEDETOMIDINE/METHADONE, White	SL-W-RABA	RABACFOSADINE, White
SL-W-MEMK	MEDETOMIDINE/METHADONE/KETAMINE, White	SL-W-RANI	RANITIDINE, White
SL-O-MELO	MELOXICAM, Orange	SL-B-REMI	REMIFENTANIL, Blue
SL-V-MELO	MELOXICAM, Violet	SL-V-RIMA	RIMADYL, Violet
SL-W-MELO	MELOXICAM, White	SL-R-ROCU	ROCURONIUM, Red
SL-G-MEPI	MEPIVACAINE, Grey	SL-V-ROCU	ROCURONIUM, Violet
SL-V-MCAT	METACAM CAT, Violet	SL-S-ROMI	ROMIFIDINE, Salmon
SL-V-META	METACAM, Violet	SL-G-ROPI	ROPIVACAINE, Grey
SL-W-META	METACAM, White	SL-W-SALB	SALBUTAMOL, White
SL-V-METR	METARAMINOL, Violet	SL-W-SA09	SALINE 0.9%, White
SL-B-METH	METHADONE, Blue	SL-W-SALF	SALINE FLUSH, White
SL-S-METO	METOCLOPRAMIDE, Salmon	SL-W-SALI	SALINE, White
SL-W-METR	METRONIDAZOLE, White	SL-W-SS09	SALINE/SODIUM CHLORIDE, White
SL-O-MIDA	MIDAZOLAM, Orange	SL-W-SEDA	SEDATION, White
SL-S-MIDA	MIDAZOLAM, Salmon	SL-B-SEDC	SEDATOR/COMFORTAN, Blue
SL-W-MITO	MITOXANTRONE, White	SL-O-SHAK	SHAKE BEFORE USE, Orange
SL-R-MIVA	MIVACURIUM, Red	SL-G-SHOR	SHORT DATE, Green
SL-O-MORP	MORPHINE N P, Orange	SL-W-SODI	SODIUM BICARBONATE, White
SL-B-MORP	MORPHINE, Blue	SL-W-SO45	SODIUM CHLORIDE 0.45%, White
SL-R-MUSC	MUSCLE RELAXANT, Red	SL-W-SO09	SODIUM CHLORIDE 0.9%, White
SL-O-NSAI	N.S.A.I.D. Orange	SL-VS-SODI	SODIUM NITROPRUSSIDE, Violet Stripe
SL-WR-NSAI	N.S.A.I.D. White (Red Writing)	SL-W-SOLU	SOLUMEDRONE, White
SL-W-NSAI	N.S.A.I.D. White	SL-W-SYNO	SYNULOX, White
SL-W-NACL	NaCl 0.9%, White	SL-W-TRAN	TRANEXAMIC ACID, White
SL-BS-NALO	NALOXONE, Blue Strip	SL-W-VENO	VENOUS, White
SL-RS-NEOS	NEOSTIGMINE, Red Stripe	SL-B-VETE	VETERGESIC, Blue
SL-O-NERV	NERVE BLOCK, Orange	SL-W-VETE	VETERGESIC, White
SL-V-NORA	NORADRENALINE, Violet	SL-R-VINB	VINBLASTINE, Red
SL-V-NORE	NOREPINEPHRINE, Violet	SL-W-VINB	VINBLASTINE, White
SL-R-OFFL	OFF LICENCE, Red	SL-R-VINC	VINCRISTINE, Red
SL-V-OMEP	OMEPRAZOLE, Violet	SL-W-VINC	VINCRISTINE, White
SL-W-OMEP	OMEPRAZOLE, White	SL-W-VINO	VINORELBINE, White
SL-S-ONDA	ONDANSETRON, Salmon	SL-W-WFIN	WATER FOR INJECTION, White
SL-R-OPEN	OPENED (DATE...), Red	SL-S-XYLA	XYLAZINE, Salmon
SL-O-PARA	PARACETAMOL, Orange	SL-W-ZINA	ZINACEF, White
SL-W-PARA	PARACETAMOL, White	SL-Y-ZINA	ZINACEF, Yellow
SL-RS-PENT	PENTOBARBITONE SODIUM, Red Stripe	SL-W-ACPM	CP/MEDETOMIDINE, White
SL-B-PETH	PETHIDINE, Blue	SL-W-ENRO	ENROFLOXACIN. White
SL-G-PHEN	PHENOBARBITONE, Green	SL-W-MARB	MARBOFLOXACIN. White
SL-V-PHEN	PHENYLEPHRINE, Violet	SL-O-ROBE	ROBENACOXIB. Orange
SL-W-POTA	POTTASSIUM CHLORIDE, White	SL-R-5FU	5FU, Red

Can't find the label you need? There are more available from our supplier that aren't listed here, contact us via email on [sales@acevetssupplies.co.uk](mailto:sales@acevetssupplies.co.uk) or call us on 01491 527027 to discuss.





# ACE

Veterinary Supplies



**ACE Veterinary Supplies Limited**

**Dogmore End Lane**

**Checkendon**

**Oxfordshire**

**RG8 0TB**

**Tel: +44 1491 527027**

**Email: [sales@acevetsupplies.co.uk](mailto:sales@acevetsupplies.co.uk)**

